



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
5.12.2008

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/8/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 4/10/08 and 5/8/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

5/12/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 12, day of May, 2008.

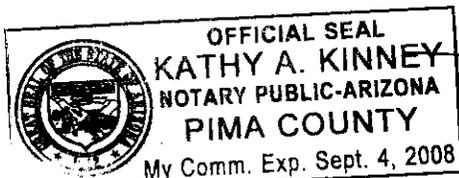

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9 SC
TBY
5.12.2008

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 5/8/08 TIME 14:47
INTOXILYZER SERIAL # 80-002803 LOCATION Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/08/2008

14:48:52

2 of 9 sc

B4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:49:36
Diagnostic Test	Pass	14:50:04
Air Blank	0.000	14:50:31
0.100 Cal Check	0.101	14:50:50
Air Blank	0.000	14:51:20
Subject Test	0.000	14:51:38
Air Blank	0.000	14:52:06
Five Minute Wait	ABT*	
Air Blank	0.000	14:53:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/08/2008

14:53:12

3 of 9 sc

BK

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:53:55
Diagnostic Test	Pass	14:54:23
Air Blank	RFI*	14:54:32
Air Blank	IPS**	14:54:40

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/08/2008

14:54:53

4 of 9 SC
PM

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:55:37
Diagnostic Test	Pass	14:56:04
Air Blank	0.000	14:56:32
0.100 Cal Check	0.100	14:56:50
Air Blank	0.000	14:57:20
Subject Test	0.087*	15:00:28
Air Blank	0.000	15:01:01
Air Blank	0.000	15:01:28
Subject Test	XXX**	15:01:44
Air Blank	0.000	15:02:17

*Deficient Sample

**Invalid Sample

Begin new deprivation period

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/08/2008

15:04:45

5 of 9 SC

BY

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:05:12
Cal Check	0.100	15:05:31
Air Blank	0.000	15:06:00

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY
5.12.2008

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/8/08 TIME 14:47
Intoxilyzer Serial # 80-002803 Location Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/08/2008 14:47:26

7 of 9 SC
D4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:47:53
Cal Check	0.101	14:48:12
Air Blank	0.000	14:48:42

Sarah Cione

Operator's Signature

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 SC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
5-12-08

Date: 5/8/08
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-002803
 New Lot Number: 735401G
 New Expiration Date: 12/28/2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/08/2008 14:45:44


9 of 9 SC

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

B4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:46:12
Cal Check	0.101	14:46:31
Air Blank	0.000	14:47:00

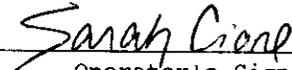

Operator's Signature

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.