



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
12-24-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-18-08 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 11-20-08 and 12-18-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


12-23-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of December, 2008.

A.M. Hernandez

Notary Public

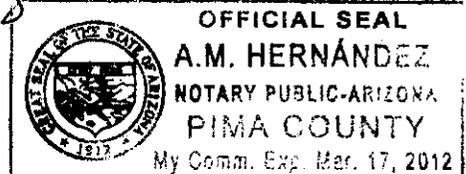

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1007  
an

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By  
12-24-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-18-08 TIME \_\_\_\_\_

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO \_\_\_\_\_

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
12/18/2008

08:32:49

2567  
an

Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517

TSY

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:33:39
Diagnostic Test	Pass	08:34:06
Air Blank	0.000	08:34:34
0.100 Cal Check	0.100	08:34:53
Air Blank	0.000	08:35:22
Subject Test	0.000	08:35:38
Air Blank	ABT*	08:35:44
Air Blank	0.000	08:36:13

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
12/18/2008

08:36:18

347  
Cu  
TS4

Standard Lot#: 804302I  
Last Changed By: C.- CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:36:47
Air Blank	IPS**	08:36:58

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
12/18/2008

08:37:04

4567  
Cu  
B4

Standard Lot#: 804302I  
Last Changed By: C.- CONLEY #34517

OAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SOAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:37:53
Diagnostic Test	Pass	08:38:21
Air Blank	0.000	08:38:48
0.100 Cal Check	0.100	08:39:07
Air Blank	0.000	08:39:36
Subject Test	0.092*	08:42:45
Air Blank	0.000	08:43:17
Air Blank	0.000	08:43:45
Subject Test	XXX**	08:44:00
Air Blank	0.000	08:44:34

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
12/18/2008

08:45:02

-----  
Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:45:29
Cal Check	0.100	08:45:48
Air Blank	0.000	08:46:17

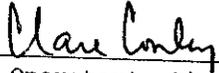  
-----  
Operator's Signature

5067  
Cu  
TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

647  
Cu

TSY  
12-24-08

QA Specialist Clare Conley AGENCY Tucson Police Department

DATE 12-18-08 TIME 0830

Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
12/18/2008 08:30:42

-----  
Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:31:10
Cal Check	0.100	08:31:28
Air Blank	0.000	08:31:58

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

7867  
CW  
TS4