



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5-11-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-7-09 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 4-9-09 and 5-7-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


5-8-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8th, day of May, 2009.

A.M. Hernández

Notary Public

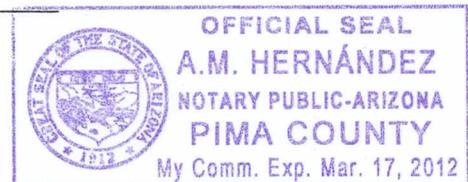

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567
Cu
Bh
5.11.09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 5-7-09 TIME 0807
INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287

cu

BM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/07/2009

08:09:30

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:10:21
Diagnostic Test	Pass	08:10:48
Air Blank	0.000	08:11:16
0.100 Cal Check	0.099	08:11:34
Air Blank	0.000	08:12:04
Subject Test	0.000	08:12:20
Air Blank	ABT*	08:12:27
Air Blank	RFI**	08:12:35

*Sequence Aborted
**RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/07/2009

08:12:43

357

cu

ASh

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	08:13:19
Air Blank	0.000	08:13:47

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/07/2009 08:13:52

4567
Ou

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517
QAS: C. CONLEY #34517
TUCSON PD
Operator: C. CONLEY #34517
TUCSON PD

BH

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:14:44
Diagnostic Test	Pass	08:15:11
Air Blank	0.000	08:15:38
0.100 Cal Check	0.099	08:15:57
Air Blank	0.000	08:16:27
Subject Test	0.065*	08:19:35
Air Blank	0.000	08:20:08
Air Blank	0.000	08:20:35
Subject Test	XXX**	08:21:05
Air Blank	0.000	08:21:39

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/07/2009 08:21:52

Standard Lot#: 8043021
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:22:19
Cal Check	0.098	08:22:38
Air Blank	0.000	08:23:08

Clare Conley

Operator's Signature

5 of 7
CH
BM

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6567

CU

BH 5.11.09

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 5-7-09 TIME 0807
Intoxilyzer Serial # 80-002503 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/07/2009

08:07:27

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:07:55
Cal Check	0.099	08:08:14
Air Blank	0.000	08:08:43

Clare Conley

Operator's Signature

7567
cu
BM