



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
6.9.2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80 - 002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/4/09 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 5/7/09 and 6/4/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/9/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of June, 20 09.

A.M. Hernandez

Notary Public

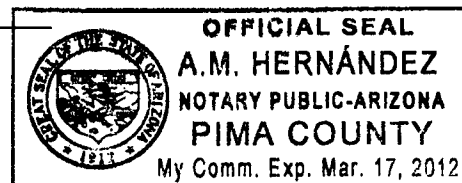

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
6-8-2009QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 6/4/09 TIME 10:03INTOXILYZER SERIAL # 80-002803 LOCATION Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

---

---

---

---

SIGNATURE

Sarah Cione

2 of 7sc

INTOXILYZER 8000

Location: ARIZONA

Serial Number : 80-002803

06/04/2009

10:04:43

154

Standard Lot#: 804302I

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:05:26
Diagnostic Test	Pass	10:05:53
Air Blank	0.000	10:06:21
0.100 Cal Check	0.099	10:06:40
Air Blank	0.000	10:07:09
Subject Test	0.000	10:07:29
Air Blank	0.000	10:07:57
Five Minute Wait	ABT*	
Air Blank	0.000	10:08:31

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
06/04/2009

10:08:36

64

Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:09:09
Air Blank	IPS**	10:09:18

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>  
104

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
06/04/2009

10:09:22

Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:10:04
Diagnostic Test	Pass	10:10:32
Air Blank	0.000	10:10:59
0.100 Cal Check	0.099	10:11:18
Air Blank	0.000	10:11:47
Subject Test	0.441*	10:14:56
Air Blank	0.000	10:15:31
Air Blank	0.000	10:15:59
Subject Test	XXX**	10:16:14
Air Blank	0.000	10:16:48

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7<sub>sc</sub>

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

06/04/2009

10:17:04

-----  
Standard Lot#: 804302I

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:17:31
Cal Check	0.098	10:17:50
Air Blank	0.000	10:18:19

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7<sub>sc</sub>

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

154  
6.8.2009

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/4/09 TIME 10:03  
Intoxilyzer Serial # 80-002803 Location Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
06/04/2009

10:03:21

-----  
Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:03:48
Cal Check	0.099	10:04:07
Air Blank	0.000	10:04:37

*Sarah Cione*

Operator's Signature