



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
5-4-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number SD-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-29-10 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 4-1-10 and 4-29-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4-30-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup>, day of April, 20 10.

A.M. Hernandez

Notary Public

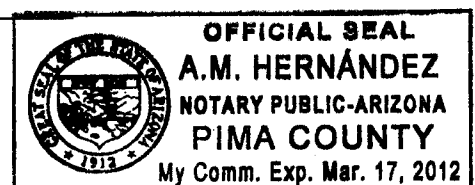

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4-29-10 TIME 0838

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
04/29/2010

08:40:03

Standard Lot#: 826001I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:51
Diagnostic Test	Pass	08:41:18
Air Blank	0.000	08:41:46
0.100 Cal Check	0.097	08:42:05
Air Blank	0.000	08:42:34
Subject Test	0.000	08:42:58
Air Blank	ABT*	08:43:04
Air Blank	0.000	08:43:33

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
04/29/2010

08:43:37

3587  
Cu  
154

Standard Lot#: 826001I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:44:24
Diagnostic Test	Pass	08:44:52
Air Blank	RFI*	08:45:02
Air Blank	IPS**	08:45:11

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
04/29/2010

08:45:17

Standard Lot#: 826001I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:46:06
Diagnostic Test	Pass	08:46:33
Air Blank	0.000	08:47:01
0.100 Cal Check	0.097	08:47:19
Air Blank	0.000	08:47:49
Subject Test	0.246*	08:50:58
Air Blank	0.000	08:51:31
Air Blank	0.000	08:51:58
Subject Test	XXX**	08:52:27
Air Blank	0.000	08:53:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

04/29/2010

08:53:21

-----  
Standard Lot#: 826001I

Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:53:48
Cal Check	0.097	08:54:07
Air Blank	0.000	08:54:37

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

657  
Ch  
TS4  
5-4-2010

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 4-29-10 TIME 0838  
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
04/29/2010

08:38:42

-----  
Standard Lot#: 826001I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:39:09
Cal Check	0.097	08:39:28
Air Blank	0.000	08:39:58

*C. Conley*

\_\_\_\_\_  
Operator's Signature