



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
6-25-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-24-10 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 5-26-10 and 6-24-10
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A

Clare Conley

Clare L. Conley


6-25-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th day of June, 2010.

Rosa Reyna

Notary Public

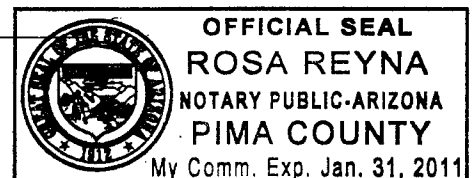

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567
en
TSY
6-25-2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 6-24-10 TIME 0838

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
06/24/2010

08:40:05

Standard Lot#: 826001I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:50
Diagnostic Test	Pass	08:41:17
Air Blank	0.000	08:41:45
0.100 Cal Check	0.098	08:42:03
Air Blank	0.000	08:42:33
Subject Test	0.000	08:42:59
Air Blank	0.000	08:43:26
Five Minute Wait	ABT*	
Air Blank	0.000	08:44:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
06/24/2010

08:44:14

Standard Lot#: 826001I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:45:33
Diagnostic Test	Pass	08:46:00
Air Blank	0.000	08:46:28
0.100 Cal Check	0.097	08:46:47
Air Blank	0.000	08:47:16
Subject Test	RFI*	08:47:26
Air Blank	IPS**	08:47:35

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
06/24/2010

08:47:44

Standard Lot#: 826001I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SOAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:48:29
Diagnostic Test	Pass	08:48:56
Air Blank	0.000	08:49:24
0.100 Cal Check	0.098	08:49:43
Air Blank	0.000	08:50:12
Subject Test	0.150*	08:53:21
Air Blank	0.000	08:53:54
Air Blank	0.000	08:54:21
Subject Test	XXX**	08:54:54
Air Blank	0.000	08:55:28

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

06/24/2010

08:55:49

Standard Lot#: 826001I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:56:17
Cal Check	0.097	08:56:35
Air Blank	0.000	08:57:05

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-24-10 TIME 0838
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

06/24/2010

08:38:32

75670a

Standard Lot#: 826001I

Last Changed By: T. GALLEGOS #28847

154

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:38:59
Cal Check	0.098	08:39:18
Air Blank	0.000	08:39:48

Clare Conley

Operator's Signature