



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
5.31.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number

80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.096 AC Cal 2 0.097 AC

Between the dates of 4/28/11 and 5/26/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -


Sarah Cione

Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30 day of

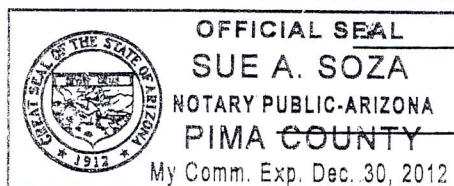

MAY, 20 11

Sue A Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
5-31-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/26/11 TIME 1346

INTOXILYZER SERIAL # 80-002803 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/26/2011

13:48:05

TS4

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:48:51
Diagnostic Test	Pass	13:49:19
Air Blank	0.000	13:49:46
0.100 Cal Check	0.098	13:50:05
Air Blank	0.000	13:50:34
Subject Test	0.000	13:50:54
Air Blank	0.000	13:51:22
Five Minute Wait	ABT*	
Air Blank	0.000	13:52:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/26/2011

13:52:12

3 of 7_{sc}

B4

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:52:43
Air Blank	IPS**	13:52:51

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
05/26/2011

13:53:47

4 of 7_{sc}
Bh

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:54:29
Diagnostic Test	Pass	13:54:57
Air Blank	0.000	13:55:24
0.100 Cal Check	0.097	13:55:43
Air Blank	0.000	13:56:12
Subject Test	0.026*	13:59:21
Air Blank	0.000	13:59:53
Air Blank	0.000	14:00:20
Subject Test	XXX**	14:00:34
Air Blank	0.000	14:01:08

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

05/26/2011

14:04:54

Standard Lot#: 719919

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:05:21
Cal Check	0.097	14:05:40
Air Blank	0.000	14:06:10

Sarah Cione

Operator's Signature

5 of 7_{sc}

BS4

EXHIBIT G-2

6 of 75c

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

124
5.31.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/26/11 TIME 1346
Intoxilyzer Serial # 80-002803 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
05/26/2011

13:46:33

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:47:00
Cal Check	0.096	13:47:19
Air Blank	0.000	13:47:48

Sarah Cione
Operator's Signature

7 of 7_{sc}

184