



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
6-27-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/9/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 5/26/11 and 6/9/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 6/9/11  
The instrument was returned to service on \_\_\_\_\_

Sarah Cione

Sarah L. Cione

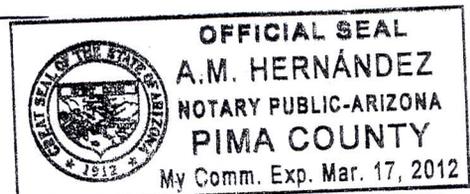
6/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27 day of June, 20 11.

A.M. Hernandez

Notary Public



\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
6-27-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 6/9/11 TIME 1535  
INTOXILYZER SERIAL # 80-002803 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
06/09/2011

15:36:30

2 of 9 sc

SM

---

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:37:16
Diagnostic Test	Pass	15:37:44
Air Blank	0.000	15:38:11
0.100 Cal Check	0.097	15:38:30
Air Blank	0.000	15:39:00
Subject Test	0.000	15:39:21
Air Blank	0.000	15:39:48
Five Minute Wait	ABT*	
Air Blank	0.000	15:40:42

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
06/09/2011

15:40:47

3 of 9 sc

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:41:32
Diagnostic Test	Pass	15:42:00
Air Blank	0.000	15:42:27
0.100 Cal Check	RFI*	15:42:40
Air Blank	IPS**	15:42:49

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
06/09/2011

15:43:11

4 of 9sc

134

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:43:54
Diagnostic Test	Pass	15:44:22
Air Blank	0.000	15:44:49
0.100 Cal Check	0.097	15:45:08
Air Blank	0.000	15:45:38
Subject Test	0.070*	15:48:46
Air Blank	0.000	15:49:19
Air Blank	0.000	15:49:46
Subject Test	XXX**	15:50:22
Air Blank	0.000	15:50:57

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
06/09/2011 15:52:47

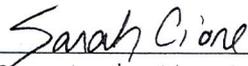
5 of 9 SC

134

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:53:14
Cal Check	0.097	15:53:33
Air Blank	0.000	15:54:02

  
\_\_\_\_\_  
Operator's Signature

6 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
6-27-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/9/11 TIME 1535  
Intoxilyzer Serial # 80-002803 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
06/09/2011 15:35:02

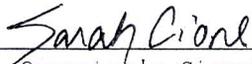
7 of 9 SC

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

TS4

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:35:30
Cal Check	0.097	15:35:48
Air Blank	0.000	15:36:18

  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

TJG  
6-27-2011

Date: 6/9/11

QAS: Sarah Ciont

Intoxilyzer Serial Number: 80-002803

Reason For Service: "Interferent Detect" occurs during subject tests when alcohol is present.

Date Out of Service: 6/9/11

Calibration Check completed prior to removal from service

Date: 6/9/11

SQAP completed prior to removal from service

Date: 6/9/11

Corrective Action:  Sent to CMI, Inc.

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_



# Service Evaluation Form

9 of 9

This form **MUST** be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.

154

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)  
 Name Terry Gallegos Phone: (520) 837-7677  
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:  
Finance / Accounts Payable Tucson Police Department Crime Lab  
Po Box 27450 270 S. Stone Ave.  
Tucson, AZ 85726-7450 Tucson, AZ 85701-1917

3. Serial Number: 80-002803 Instrument Model: I8000

4. Detailed Description of Problem:  
When subject test performed (+ alcohol is present), interferent detect occurs.  
Completes alcohol free test fine.

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\*

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator  
 Name (Please Print) Title  
[Signature] 6.27.2011  
 Signature Date

Ship to:

CMI, Inc.  
 Attn: Service Dept.  
 316 East Ninth Street  
 Owensboro, KY 42303

No, please send estimate before repairs are made.  
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.