



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*By  
2.23.2012*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/7/12 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of - and - the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on 2/7/12.

Sarah Cione

Sarah L. Cione


2/10/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup>, day of February, 20 12.

Elma Escarcega

Notary Public

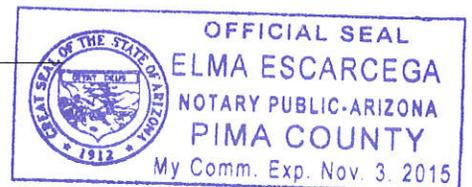

\*Includes documents specified in A.A.C. R13-10-104.A.5



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 15<sub>sc</sub>

T54  
2.23.2012

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002803

Date Returned to Service: 2/7/12

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSh  
2-23-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/7/12 TIME 1435

INTOXILYZER SERIAL # 80-002803 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS Returned from service in Sept, 2011. Had to go back to CMI because it would not pass initial diagnostics

SIGNATURE Sarah Cione

3 of 15 sc

Tsh

INTOXILYZER 8000

Location:

Serial Number: 80-002803

Core Version: 8105.48

02/07/2012

14:37:25

---

Standard Lot#: 23511100A2

Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616

TUCSON PD

---

Operator: S. CIONE #50616

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:38:08
Diagnostic Test	Pass	14:38:43
Air Blank	0.000	14:39:10
0.100 Cal Check	0.099	14:39:31
Air Blank	0.000	14:40:01
Subject Test	0.000	14:40:14
Air Blank	0.000	14:40:41
Wait	ABT*	
Air Blank	0.000	14:41:16

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

4 of 15<sub>sc</sub>

INTOXILYZER 8000

Location:

Serial Number: 80-002803

Core Version: 8105.48

02/07/2012

14:41:22

T34

---

Standard Lot#: 23511100A2

Last Changed By: S. CIONE

#50616

---

QAS: S. CIONE

#50616

TUCSON PD

---

Operator: S. CIONE

#50616

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	14:41:51
Air Blank	IPS**	14:41:58

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

5 of 15 sc

INTOXILYZER 8000

Location:

Serial Number: 80-002803

Core Version: 8105.48

02/07/2012

14:42:04

TS4

Standard Lot#: 23511100A2

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:42:47
Diagnostic Test	Pass	14:43:22
Air Blank	0.000	14:43:49
0.100 Cal Check	0.100	14:44:10
Air Blank	0.000	14:44:40
Subject Test	0.076*	14:47:47
Air Blank	0.000	14:48:20
Air Blank	0.000	14:48:48
Subject Test	XXX**	14:49:02
Air Blank	0.000	14:49:36

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

6 of 15 sc

Sh

INTOXILYZER 8000  
Location:  
Serial Number: 80-002803  
Core Version: 8105.48  
02/07/2012 14:50:20

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:50:49
0.100 Cal Check	0.099	14:51:09
Air Blank	0.000	14:51:39

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

7 of 15<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4  
2.23.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/7/12 TIME 1435  
Intoxilyzer Serial # 80-002803 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

8 of 15<sub>sc</sub>

TSh

INTOXILYZER 8000  
Location:  
Serial Number: 80-002803  
Core Version: 8105.48  
02/07/2012 14:35:51

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:36:20
0.100 Cal Check	0.100	14:36:41
Air Blank	0.000	14:37:11

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 15<sub>sc</sub>

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TS4  
2.23.2012

Date: 2/7/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002803

New Lot Number: 23511100A2

New Expiration Date: 11/1/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

10 of 15<sub>sc</sub>  
TSH

INTOXILYZER 8000  
Location:  
Serial Number: 80-002803  
Core Version: 8105.48  
02/07/2012 14:34:23

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:34:51
0.100 Cal Check	0.100	14:35:12
Air Blank	0.000	14:35:42

Preventative Maintenance  
Performed -- Timer Reset



12 of 15 sc

T24

# Certificate of Calibration

## **Intoxilyzer**® 8000

*This is to certify that the calibration of **Intoxilyzer**® serial number 80-002803, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.*

Date Sept. 21, 2011 Signed *Debbie Palmer*  
Technician



316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



# Service Evaluation Form

13 of 15 SC

Ph

This form **MUST** be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)  
 Name Terry Gallegos Phone: (520) 837-7677  
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:  
Finance / Accounts Payable Tucson Police Department Crime Lab  
Po Box 27450 1306 W. Miracle Mile  
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-002803 Instrument Model: I8000

4. Detailed Description of Problem:  
Came back from service. Upon 1st time plugging in, will not pass Diagnostics -  
Analytical Stability Fails.

**\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\***

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator  
 Name (Please Print) Title  
Terry Gallegos 11.22.2011  
 Signature Date

Ship to:  
 CMI, Inc.  
 Attn: Service Dept.  
 316 East Ninth Street  
 Owensboro, KY 42303

No, please send estimate before repairs are made.  
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.



15 of 15 sc

B4

# Certificate of Calibration

## **Intoxilyzer**® 8000

This is to certify that the calibration of **Intoxilyzer**® serial number 80-002803, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

Date Jan. 27, 2012 Signed Dale Rohrer

Technician

**CMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.