



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
6-27-2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/21/12 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 5/24/12 and 6/21/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/21/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21<sup>st</sup>, day of June, 2012.

Elma Escarcega

Notary Public

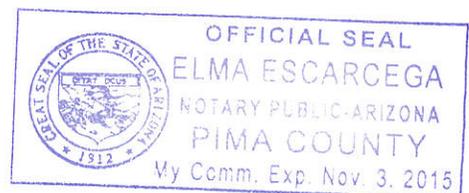

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
6-27-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/21/12 TIME 949

INTOXILYZER SERIAL # 80-002803 LOCATION ODW-Tm 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7<sub>sc</sub>

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/21/2012

09:51:01

TS4

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

OAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:51:45
Diagnostic Test	Pass	09:52:20
Air Blank	0.000	09:52:47
0.100 Cal Check	0.098	09:53:08
Air Blank	0.000	09:53:38
Subject Test	0.000	09:53:54
Air Blank	0.000	09:54:21
Wait	ABT*	
Air Blank	0.000	09:54:55

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/21/2012 09:55:01

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:55:36
Air Blank	IPS**	09:55:44

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 sc

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/21/2012

09:55:53

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:56:37
Diagnostic Test	Pass	09:57:12
Air Blank	0.000	09:57:39
0.100 Cal Check	0.098	09:58:00
Air Blank	0.000	09:58:29
Subject Test	0.068*	10:01:37
Air Blank	0.000	10:02:09
Air Blank	0.000	10:02:36
Subject Test	XXX**	10:03:04
Air Blank	0.000	10:03:38

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7<sub>sc</sub>

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/21/2012 10:03:47

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:04:16
0.100 Cal Check	0.098	10:04:36
Air Blank	0.000	10:05:06

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
6-27-2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/21/12 TIME 949  
Intoxilyzer Serial # 80-002803 Location ODW-Tm 2

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7sc

154

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/21/2012 09:49:27

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:49:56
0.100 Cal Check	0.099	09:50:16
Air Blank	0.000	09:50:46

Preventative Maintenance  
Performed -- Timer Reset