



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TSh  
 9.20.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/14/12 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 8/16/12 and 9/14/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/17/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of September, 2012.

Elma Escarcega

Notary Public

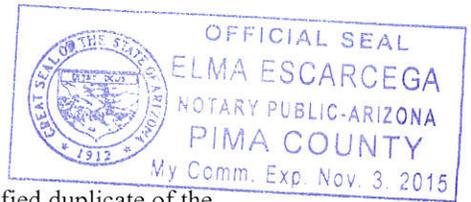

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

TS4  
9.20.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 9/14/12 TIME 1005  
INTOXILYZER SERIAL # 80-002803 LOCATION Tm 2 - ODW

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7 sc

TS4

INTOXILYZER 8000  
 Location: TUCSON PD  
 Serial Number: 80-002803  
 Core Version: 8105.48  
 09/14/2012 10:07:21

---

Standard Lot#: 23511100A2  
 Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
 TUCSON PD

---

Operator: S. CIONE #50616  
 TUCSON PD

---

Subject: SQAP, SQAP  
 DOB: 11/11/1980  
 Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:08:05
Diagnostic Test	Pass	10:08:40
Air Blank	0.000	10:09:07
0.100 Cal Check	0.097	10:09:28
Air Blank	0.000	10:09:57
Subject Test	0.000	10:10:14
Air Blank	0.000	10:10:41
Wait	ABT*	
Air Blank	0.000	10:11:18

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>sc</sub>

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/14/2012 10:11:23

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:12:07
Diagnostic Test	Pass	10:12:42
Air Blank	0.000	10:13:09
0.100 Cal Check	0.097	10:13:30
Air Blank	0.000	10:14:00
Subject Test	RFI*	10:14:13
Air Blank	IPS**	10:14:22

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/14/2012 10:14:28

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:15:09
Diagnostic Test	Pass	10:15:44
Air Blank	0.000	10:16:11
0.100 Cal Check	0.098	10:16:32
Air Blank	0.000	10:17:02
Subject Test	0.139*	10:20:09
Air Blank	0.000	10:20:41
Air Blank	0.000	10:21:09
Subject Test	XXX**	10:21:25
Air Blank	0.000	10:22:00

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/14/2012 10:22:11

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:22:39
0.100 Cal Check	0.096	10:23:00
Air Blank	0.000	10:23:30

Preventative Maintenance  
Performed -- Timer Reset

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

TS4

9.20.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/14/12 TIME 1005  
Intoxilyzer Serial # 80-002803 Location Tm2-ODW

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/14/2012 10:05:47

7 of 7  
sc  
TS4

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:06:15
0.100 Cal Check	0.096	10:06:36
Air Blank	0.000	10:07:05

Preventative Maintenance  
Performed -- Timer Reset