



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSh
7-24-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/18/13 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 6/21/13 and 7/18/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/22/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22nd, day of July, 2013.

Elma Escarcega

Notary Public

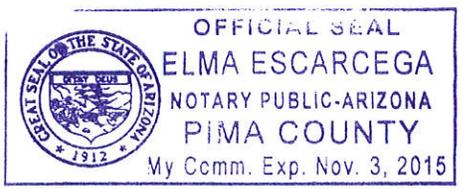

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
7.24.2013

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/18/13 TIME 935

INTOXILYZER SERIAL # 80-002803 LOCATION ODW-Tm 2

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 85c

TBY

INTOXILYZER 8000
 Location: TUCSON PD
 Serial Number: 80-002803
 Core Version: 8105.48
 07/18/2013 09:37:30

Standard Lot#: 30412100A1
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:38:13
Diagnostic Test	Pass	09:38:48
Air Blank	0.000	09:39:16
0.100 Cal Check	0.097	09:39:36
Air Blank	0.000	09:40:06
Subject Test	0.000	09:40:36
Air Blank	0.000	09:41:03
Wait	ABT*	
Air Blank	0.000	09:41:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8_{SC}

TD4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
07/18/2013 09:41:49

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:42:32
Diagnostic Test	Pass	09:43:07
Air Blank	0.000	09:43:35
0.100 Cal Check	RFI*	09:43:51
Air Blank	IPS**	09:43:59

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
07/18/2013

09:44:04

4 of 8 sc
By

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:44:46
Diagnostic Test	Pass	09:45:21
Air Blank	0.000	09:45:48
0.100 Cal Check	0.097	09:46:09
Air Blank	0.000	09:46:38
Subject Test	RFI*	09:49:14
Air Blank	0.000	09:49:47

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
07/18/2013

09:50:18

5 of 8 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:51:01
Diagnostic Test	Pass	09:51:36
Air Blank	0.000	09:52:04
0.100 Cal Check	0.096	09:52:24
Air Blank	0.000	09:52:54
Subject Test	0.185*	09:56:01
Air Blank	0.000	09:56:33
Air Blank	0.000	09:57:00
Subject Test	XXX**	09:57:24
Air Blank	0.000	09:57:58

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
07/18/2013 09:58:08

6 of 8_{SC}
TBH

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:58:36
0.100 Cal Check	0.096	09:58:57
Air Blank	0.000	09:59:27

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8
JC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSJ
7.24.2013

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/18/13 TIME 935
Intoxilyzer Serial # 80-002803 Location ODW-Tm 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
07/18/2013 09:35:50

8 of 8
JC
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:36:19
0.100 Cal Check	0.097	09:36:40
Air Blank	0.000	09:37:09

Preventative Maintenance
Performed -- Timer Reset