



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
6.19.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/19/14 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 5/22/14 and 6/19/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/19/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th, day of June, 2014.

Elma Escarcega

Notary Public

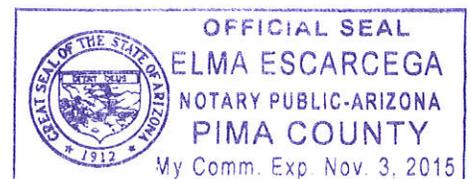

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 8 sc  
TS4  
6-19-2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/19/14 TIME 756

INTOXILYZER SERIAL # 80-002803 LOCATION ODW-Tm2

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 8 sc  
TSh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

07:57:48

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:58:31
Diagnostic Test	Pass	07:59:06
Air Blank	0.000	07:59:33
0.100 Cal Check	0.098	07:59:54
Air Blank	0.000	08:00:23
Subject Test	0.000	08:00:45
Air Blank	0.000	08:01:13
Wait	ABT*	
Air Blank	0.000	08:01:48

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

08:01:54

3 of 8 sc  
B4

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:02:40
Diagnostic Test	Pass	08:03:15
Air Blank	0.000	08:03:42
0.100 Cal Check	RFI*	08:03:57
Air Blank	IPS**	08:04:07

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

08:04:13

4 of 8 sc  
TJ4

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:05:17
Diagnostic Test	Pass	08:05:52
Air Blank	0.000	08:06:20
0.100 Cal Check	0.097	08:06:41
Air Blank	0.000	08:07:10
Subject Test	RFI*	08:09:12
Air Blank	0.000	08:09:44

\*RFI Detect

- unknown source of RFI

**Not a Successfully Completed Test Sequence**

5 of 8 sc  
TBH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

08:10:26

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:11:12
Diagnostic Test	Pass	08:11:47
Air Blank	0.000	08:12:14
0.100 Cal Check	0.097	08:12:35
Air Blank	0.000	08:13:04
Subject Test	0.179*	08:16:12
Air Blank	0.000	08:16:44
Air Blank	0.000	08:17:11
Subject Test	XXX**	08:17:27
Air Blank	0.000	08:18:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

08:18:17

6 of 8<sub>sc</sub>  
By

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:18:46
0.100 Cal Check	0.096	08:19:06
Air Blank	0.000	08:19:36

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8  
sc  
TS4  
6.19.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/19/14 TIME 756  
Intoxilyzer Serial # 80-002803 Location ODW-Tm 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/19/2014

07:56:07

8 of 8 sc  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:56:35
0.100 Cal Check	0.098	07:56:56
Air Blank	0.000	07:57:25

Preventative Maintenance  
Performed -- Timer Reset