



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS6  
8.4.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/17/14 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 6/19/14 and 7/17/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/18/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 18<sup>th</sup>, day of July, 2014.

Elma Escarcega

Notary Public

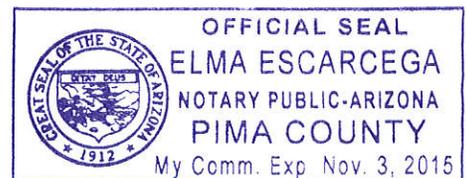

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS  
8.4.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/17/14 TIME 914

INTOXILYZER SERIAL # 80-002803 LOCATION ODW-Tm 2

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/17/2014

09:16:08

2 of 7  
B4

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:16:51 |
| Diagnostic Test | Pass   | 09:17:26 |
| Air Blank       | 0.000  | 09:17:53 |
| 0.100 Cal Check | 0.096  | 09:18:14 |
| Air Blank       | 0.000  | 09:18:43 |
| Subject Test    | 0.000  | 09:19:24 |
| Air Blank       | 0.000  | 09:19:51 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 09:20:28 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/17/2014

09:20:44

3077  
SC  
B4

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:21:24 |
| Diagnostic Test | Pass   | 09:21:59 |
| Air Blank       | 0.000  | 09:22:27 |
| 0.100 Cal Check | 0.096  | 09:22:47 |
| Air Blank       | 0.000  | 09:23:17 |
| Subject Test    | RFI*   | 09:23:29 |
| Air Blank       | IPS**  | 09:23:37 |

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/17/2014

09:23:44

4 of 7  
SC  
B4

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:24:26 |
| Diagnostic Test | Pass   | 09:25:01 |
| Air Blank       | 0.000  | 09:25:28 |
| 0.100 Cal Check | 0.096  | 09:25:49 |
| Air Blank       | 0.000  | 09:26:19 |
| Subject Test    | 0.227* | 09:29:26 |
| Air Blank       | 0.000  | 09:29:58 |
| Air Blank       | 0.000  | 09:30:26 |
| Subject Test    | XXX**  | 09:30:45 |
| Air Blank       | 0.000  | 09:31:18 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/17/2014 09:31:48

5 of 7  
sc  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:32:16 |
| 0.100 Cal Check | 0.095  | 09:32:37 |
| Air Blank       | 0.000  | 09:33:07 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

607 7 SC  
B4  
8.4.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 7/17/14 TIME 914

Intoxilyzer Serial # 80-002803 Location ODW - Tm 2

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/17/2014 09:14:19

7 of 7 sc  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

-----  
QAS: S. CIONE #50616  
TUCSON PD

-----  
Operator: S. CIONE #50616  
TUCSON PD

-----  
Test g/210L Time  
-----  
Air Blank 0.000 09:14:47  
0.100 Cal Check 0.096 09:15:08  
Air Blank 0.000 09:15:38

Preventative Maintenance  
Performed -- Timer Reset