



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3.12.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/5/2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 2/4/2015 and 3/5/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis

3/11/2015


ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 11th, day of March, 2015.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

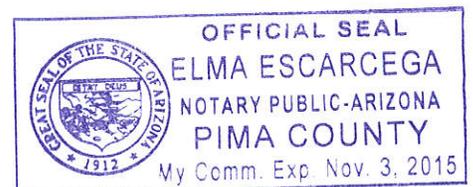

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
Bh
3.12.2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 3.5.2015 TIME 14:33

INTOXILYZER SERIAL # 80-002803 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/05/2015

14:35:18

2 of 7
WLL
B4

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:36:06
Diagnostic Test	Pass	14:36:41
Air Blank	0.000	14:37:09
0.100 Cal Check	0.100	14:37:29
Air Blank	0.000	14:37:59
Subject Test	0.000	14:38:26
Air Blank	0.000	14:38:54
Wait	ABT*	
Air Blank	0.000	14:39:32

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/05/2015

14:40:09

3 of 7
WV
TSH

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:41:19
Diagnostic Test	Pass	14:41:54
Air Blank	0.000	14:42:21
0.100 Cal Check	RFI*	14:42:38
Air Blank	IPS**	14:42:48

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/05/2015

14:43:07

4 of 7
WLL
Bh

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:44:17
Diagnostic Test	Pass	14:44:52
Air Blank	0.000	14:45:20
0.100 Cal Check	0.100	14:45:41
Air Blank	0.000	14:46:10
Subject Test	0.438*	14:49:18
Air Blank	0.000	14:49:59
Air Blank	0.000	14:50:26
Subject Test	XXX**	14:50:47
Air Blank	0.000	14:51:21

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/05/2015 14:51:50

5 of 7
WLL
BS

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:52:20
0.100 Cal Check	0.098	14:52:41
Air Blank	0.000	14:53:12

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TJH
3.12.2015

QA Specialist WILL LOOMIS 53545 AGENCY Tucson Police Department
DATE 3.5.2015 TIME 14:33
Intoxilyzer Serial # 80-002803 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/05/2015 14:33:10

7 of 7
TS4
WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:33:39
0.100 Cal Check	0.100	14:34:00
Air Blank	0.000	14:34:30

Preventative Maintenance
Performed -- Timer Reset