



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

Handwritten note: B4 4.10.2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/1/2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 3/5/2015 and 4/1/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on [blank]. The instrument was returned to service on [blank].

Handwritten signature: William L Loomis

William L. Loomis


4/2/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2nd day of April, 2015.

Handwritten signature: Elma Escarcega

Notary Public

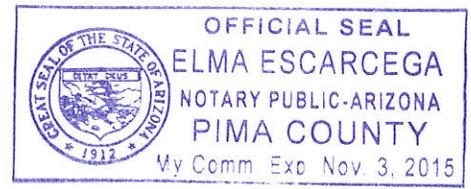

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

*1 of 7 will
Bq 4.10.2015*

QA SPECIALIST Will Loomis 53545 AGENCY Tucson Police Department

DATE 4/1/2015 TIME 16:40

INTOXILYZER SERIAL # 80-002803 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
04/01/2015

16:42:31

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:43:23
Diagnostic Test	Pass	16:43:57
Air Blank	0.000	16:44:25
0.100 Cal Check	0.100	16:44:46
Air Blank	0.000	16:45:15
Subject Test	0.000	16:45:35
Air Blank	0.000	16:46:02
Wait	ABT*	
Air Blank	0.000	16:46:38

*Sequence Aborted

2 of 7
wrc
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
04/01/2015

16:46:55

3 of 7
WLL
TSH

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:47:51
Diagnostic Test	Pass	16:48:26
Air Blank	0.000	16:48:53
0.100 Cal Check	RFI*	16:49:07
Air Blank	IPS**	16:49:16

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
04/01/2015

16:49:23

4 of 7
TSG WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:50:14
Diagnostic Test	Pass	16:50:49
Air Blank	0.000	16:51:16
0.100 Cal Check	0.100	16:51:37
Air Blank	0.000	16:52:07
Subject Test	0.173*	16:55:14
Air Blank	0.000	16:55:48
Air Blank	0.000	16:56:16
Subject Test	XXX**	16:56:39
Air Blank	0.000	16:57:13

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
04/01/2015 16:57:38

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:58:13
0.100 Cal Check	0.098	16:58:34
Air Blank	0.000	16:59:03

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLL
B4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 4/1/2015 TIME 16:40

Intoxilyzer Serial # 80-002803 Location CRIME LAB

6 of 7 WLL
TS4
4.10.2015

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
04/01/2015 16:40:30

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:41:00
0.100 Cal Check	0.100	16:41:21
Air Blank	0.000	16:41:51

Preventative Maintenance
Performed -- Timer Reset

7 of 7
0 WLL
TS4