



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
8-6-2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/23/2015 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 6/25/2015 and 7/23/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
 The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


7/27/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 27<sup>th</sup>, day of July, 2015.

Elma Escarcega

Notary Public

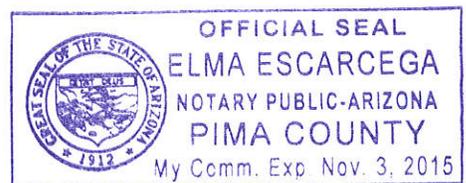

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
T04 8.6.2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 7/23/2015 TIME 14:39

INTOXILYZER SERIAL # 80-002803 LOCATION CRIME LAB

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/23/2015

14:42:06

2 of 7  
WLL  
TS4

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:42:54
Diagnostic Test	Pass	14:43:29
Air Blank	0.000	14:43:57
0.100 Cal Check	0.097	14:44:17
Air Blank	0.000	14:44:47
Subject Test	0.000	14:45:04
Air Blank	0.000	14:45:31
Wait	ABT*	
Air Blank	0.000	14:46:15

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/23/2015

14:46:43

3 of 7  
WV  
TJ

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:49:45
Diagnostic Test	Pass	14:50:20
Air Blank	0.000	14:50:47
0.100 Cal Check	RFI*	14:51:02
Air Blank	IPS**	14:51:09

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/23/2015

14:51:14

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:52:12
Diagnostic Test	Pass	14:52:47
Air Blank	0.000	14:53:14
0.100 Cal Check	0.097	14:53:35
Air Blank	0.000	14:54:04
Subject Test	0.144*	14:57:12
Air Blank	0.000	14:57:44
Air Blank	0.000	14:58:11
Subject Test	XXX**	14:58:30
Air Blank	0.000	14:59:04

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
wll  
TSh

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/23/2015 14:59:26

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:59:55
0.100 Cal Check	0.096	15:00:16
Air Blank	0.000	15:00:46

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
TS4 WLL

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TS48.6.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 7/23/2015 TIME 14:39  
Intoxilyzer Serial # 80-002803 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
07/23/2015 14:39:43

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:40:12
0.100 Cal Check	0.097	14:40:33
Air Blank	0.000	14:41:02

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
ww  
TS4