



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSJ
3.23.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/3/2016 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 2/4/2016 and 3/3/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


3/7/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 7th, day of March, 2016.

Elma Escarcega

Notary Public

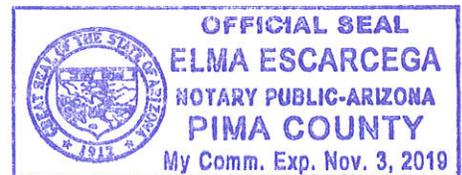

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TS4 3.23.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 3/3/2016 TIME 11:06
INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE William Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/03/2016

11:08:26

2 of 7
wlu
TS4

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:09:14
Diagnostic Test	Pass	11:09:50
Air Blank	0.000	11:10:17
0.100 Cal Check	0.097	11:10:38
Air Blank	0.000	11:11:07
Subject Test	0.000	11:11:25
Air Blank	0.000	11:11:52
Wait	ABT*	
Air Blank	0.000	11:12:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/03/2016

11:12:32

3 of 7
WLL
TSG

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:13:23
Diagnostic Test	Pass	11:13:58
Air Blank	0.000	11:14:25
0.100 Cal Check	RFI*	11:14:39
Air Blank	IPS**	11:14:48

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/03/2016

11:14:51

4 of 7
will
TS4

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:15:39
Diagnostic Test	Pass	11:16:14
Air Blank	0.000	11:16:41
0.100 Cal Check	0.097	11:17:02
Air Blank	0.000	11:17:31
Subject Test	0.261*	11:20:39
Air Blank	0.000	11:21:12
Air Blank	0.000	11:21:40
Subject Test	XXX**	11:21:55
Air Blank	0.000	11:22:30

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/03/2016 11:22:51

5 of 7
will
TS4

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:23:21
0.100 Cal Check	0.095	11:23:42
Air Blank	0.000	11:24:12

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TS4
3.23.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 3/3/2016 TIME 11:06
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/03/2016 11:06:55

7 of 7
WLL
TS4

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:07:24
0.100 Cal Check	0.097	11:07:45
Air Blank	0.000	11:08:14

Preventative Maintenance
Performed -- Timer Reset