



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy
4.7.2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/31/2016 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 3/3/2016 and 3/31/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


4/4/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2/4th, day of April, 2016.

Elma Escarcega

Notary Public

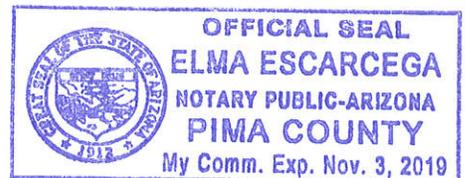

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
wn
Thu 4.7.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 3/31/2016 TIME 14:06

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Willie Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/31/2016

14:09:16

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:10:16
Diagnostic Test	Pass	14:10:51
Air Blank	0.000	14:11:18
0.100 Cal Check	0.097	14:11:39
Air Blank	0.000	14:12:09
Subject Test	0.000	14:12:25
Air Blank	0.000	14:12:53
Wait	ABT*	
Air Blank	0.000	14:13:28

*Sequence Aborted

2 of 7
WLL
TSH

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/31/2016

14:13:32

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:14:22
Diagnostic Test	Pass	14:14:57
Air Blank	0.000	14:15:24
0.100 Cal Check	RFI*	14:15:38
Air Blank	IPS**	14:15:46

*RFI Detect

**Improper Sample

3 of 7
WLL
BY

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/31/2016

14:15:58

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:16:48
Diagnostic Test	Pass	14:17:22
Air Blank	0.000	14:17:50
0.100 Cal Check	0.096	14:18:11
Air Blank	0.000	14:18:40
Subject Test	0.213*	14:21:48
Air Blank	0.000	14:22:20
Air Blank	0.000	14:22:47
Subject Test	XXX**	14:23:14
Air Blank	0.000	14:23:47

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
0 WU
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/31/2016 14:24:01

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:24:31
0.100 Cal Check	0.095	14:24:52
Air Blank	0.000	14:25:21

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLL
TSG

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLV
4.7.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 3/31/2016 TIME 14:06
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002803
Core Version: 8105.48
03/31/2016 14:06:58

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:07:27
0.100 Cal Check	0.096	14:07:48
Air Blank	0.000	14:08:17

Preventative Maintenance
Performed -- Timer Reset

7 of 7
0 WLL
TSh