



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
6-28-2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/23/2016 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 5/26/2016 and 6/23/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


6/24/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 24<sup>th</sup>, day of June, 2016.

Elma Escarcega

Notary Public

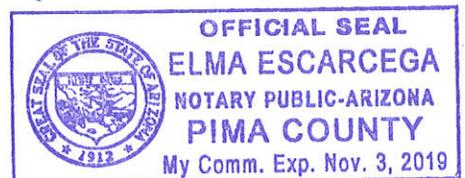

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8  
ML  
B4 6-28-2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 6/23/2016 TIME 11:12

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Willie Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016

11:14:52

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:15:44
Diagnostic Test	Pass	11:16:19
Air Blank	0.000	11:16:46
0.100 Cal Check	0.094	11:17:07
Air Blank	0.000	11:17:36
Subject Test	0.000	11:17:52
Air Blank	0.000	11:18:19
Wait	ABT*	
Air Blank	0.000	11:18:54

\*Sequence Aborted

2 of 8  
0 WLL  
TBL

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016

11:18:59

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:20:00
Diagnostic Test	Pass	11:20:35
Air Blank	0.000	11:21:03
0.100 Cal Check	RFI*	11:21:16
Air Blank	IPS**	11:21:24

\*RFI Detect  
\*\*Improper Sample

3 of 8  
0 WIL  
T84

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016

11:21:28

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:22:23
Diagnostic Test	Pass	11:22:58
Air Blank	0.000	11:23:25
0.100 Cal Check	0.095	11:23:46
Air Blank	0.000	11:24:16
Subject Test	0.240*	11:27:23
Air Blank	RFI**	11:27:39
Air Blank	0.000	11:28:08

\*Deficient Sample

\*\*RFI Detect

4 of 8  
6 WLV

TSG

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016

11:28:33

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:29:29
Diagnostic Test	Pass	11:30:04
Air Blank	0.000	11:30:31
0.100 Cal Check	0.093	11:30:52
Air Blank	0.000	11:31:22
Subject Test	0.208*	11:34:29
Air Blank	0.000	11:35:02
Air Blank	0.000	11:35:29
Subject Test	XXX**	11:35:48
Air Blank	0.000	11:36:22

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

5 of 8  
0 WLL  
TJG

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016 11:36:36

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:37:07
0.100 Cal Check	0.094	11:37:28
Air Blank	0.000	11:37:58

Preventative Maintenance  
Performed -- Timer Reset

6 of 8  
0 will  
TBA

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8  
WLL  
T84 6-28-2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 6/23/2016 TIME 11:12  
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
06/23/2016 11:12:43

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:13:13
0.100 Cal Check	0.095	11:13:34
Air Blank	0.000	11:14:03

Preventative Maintenance  
Performed -- Timer Reset

8 of 8  
0 WLL  
D4