



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*BS4  
9.22.2016*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/2016 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 8/18/2016 and 9/15/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
 The instrument was returned to service on \_\_\_\_\_.

William L Loomis


William L. Loomis  
9/19/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 19<sup>th</sup>, day of September, 2016.

Elma Escarcega

Notary Public

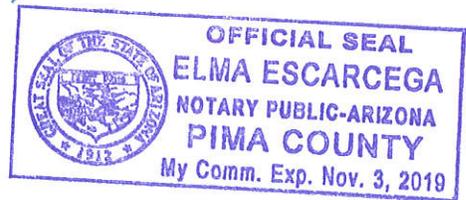

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
Bh 9.22.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 9/15/2016 TIME 11:34

INTOXILYZER SERIAL # 80-002803 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE William Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/15/2016

11:36:12

2 of 7  
WLL  
TS4

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:37:05
Diagnostic Test	Pass	11:37:40
Air Blank	0.000	11:38:07
0.100 Cal Check	0.098	11:38:28
Air Blank	0.000	11:38:58
Subject Test	0.000	11:39:12
Air Blank	0.000	11:39:39
Wait	ABT*	
Air Blank	0.000	11:40:13

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/15/2016

11:40:18

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:41:12
Diagnostic Test	Pass	11:41:47
Air Blank	0.000	11:42:14
0.100 Cal Check	RFI*	11:42:28
Air Blank	IPS**	11:42:36

\*RFI Detect

\*\*Improper Sample

3 of 7  
WLL  
TS4

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/15/2016

11:42:48

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:43:41
Diagnostic Test	Pass	11:44:16
Air Blank	0.000	11:44:43
0.100 Cal Check	0.098	11:45:04
Air Blank	0.000	11:45:33
Subject Test	0.230*	11:48:41
Air Blank	0.000	11:49:13
Air Blank	0.000	11:49:41
Subject Test	XXX**	11:49:57
Air Blank	0.000	11:50:31

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
WLL  
TSY

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/15/2016

11:50:47

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:51:16
0.100 Cal Check	0.097	11:51:37
Air Blank	0.000	11:52:06

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WLL  
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL

By 9.22.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 9/15/2016 TIME 11:34  
Intoxilyzer Serial # 80-002803 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002803  
Core Version: 8105.48  
09/15/2016 11:34:44

7 of 7  
WLL

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

BY

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:35:13
0.100 Cal Check	0.098	11:35:34
Air Blank	0.000	11:36:03

Preventative Maintenance  
Performed -- Timer Reset