



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
3/13/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006048 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-26-2015 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 1-29-2015 and 2-26-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos


3-10-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup> day of March, 2015.

Elma Escarcega

Notary Public

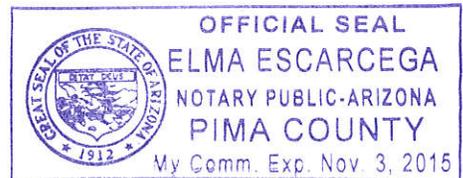

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10F7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
3/13/15

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 0528847 AGENCY Tucson Police Department

DATE 2.26.2015 TIME 1101

INTOXILYZER SERIAL # 80-006048 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

2077  
Tg  
sc

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
02/26/2015 11:03:54

---

Standard Lot#: 25414100A4  
Last Changed By: J. KNEPPER #46585

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:04:55
Diagnostic Test	Pass	11:05:30
Air Blank	0.000	11:05:57
0.100 Cal Check	0.100	11:06:18
Air Blank	0.000	11:06:48
Subject Test	0.000	11:08:22
Air Blank	0.000	11:08:50
Wait	ABT*	
Air Blank	0.000	11:09:25

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7  
TS4  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
02/26/2015 11:09:51

---

Standard Lot#: 25414100A4  
Last Changed By: J. KNEPPER #46585

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:11:35
Diagnostic Test	Pass	11:12:10
Air Blank	0.000	11:12:37
0.100 Cal Check	RFI*	11:12:54
Air Blank	IPS**	11:13:09

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7  
TJ

sc

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
02/26/2015

11:14:48

---

Standard Lot#: 25414100A4  
Last Changed By: J. KNEPPER #46585

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:15:38
Diagnostic Test	Pass	11:16:13
Air Blank	0.000	11:16:41
0.100 Cal Check	0.100	11:17:02
Air Blank	0.000	11:17:31
Subject Test	0.019*	11:20:38
Air Blank	0.000	11:21:08
Air Blank	0.000	11:21:36
Subject Test	XXX**	11:21:58
Air Blank	0.000	11:22:31

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

507  
154

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
02/26/2015 12:45:59

5c

-----  
Standard Lot#: 25414100A4  
Last Changed By: J. KNEPPER #46585  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:46:31
0.100 Cal Check	0.100	12:46:52
Air Blank	0.000	12:47:21

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6077  
TS4  
SC  
3/13/15

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos AGENCY Tucson Police Department  
DATE 2.26.2015 TIME 1101  
Intoxilyzer Serial # 80-006048 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution      AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
02/26/2015 11:01:19

7 of 7  
154

-----  
Standard Lot#: 25414100A4  
Last Changed By: J. KNEPPER #46585  
-----

SC

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:01:49
0.100 Cal Check	0.100	11:02:10
Air Blank	0.000	11:02:39

Preventative Maintenance  
Performed -- Timer Reset