



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006048 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL
 4.20.2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-26-2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2-26-2015 and 3-26-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos

4.16.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16th, day of April, 2015.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

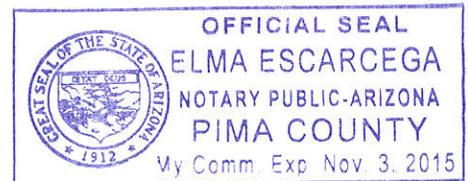


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10F7
TSh

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

WLV
4.20.2015

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 3.26.2015 TIME 1229

INTOXILYZER SERIAL # 80-006048 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
03/26/2015

12:32:25

2 of 7
TS4

WLL

Standard Lot#: 25414100A4
Last Changed By: J. KNEPPER #46585

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:33:14
Diagnostic Test	Pass	12:33:49
Air Blank	0.000	12:34:16
0.100 Cal Check	0.101	12:34:37
Air Blank	0.000	12:35:06
Subject Test	0.000	12:35:27
Air Blank	0.000	12:35:54
Wait	ABT*	
Air Blank	0.000	12:36:40

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
03/26/2015

12:36:55

3 of 7
JG

WLL

Standard Lot#: 25414100A4
Last Changed By: J. KNEPPER #46585

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:38:43
Diagnostic Test	Pass	12:39:18
Air Blank	0.000	12:39:46
0.100 Cal Check	RFI*	12:40:02
Air Blank	IPS**	12:40:15

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
03/26/2015

12:41:13

4 of 7
Dy

Standard Lot#: 25414100A4
Last Changed By: J. KNEPPER #46585

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:42:01
Diagnostic Test	Pass	12:42:36
Air Blank	0.000	12:43:03
0.100 Cal Check	0.099	12:43:24
Air Blank	0.000	12:43:53
Subject Test	0.141*	12:47:01
Air Blank	0.000	12:47:33
Air Blank	0.000	12:48:00
Subject Test	XXX**	12:48:16
Air Blank	0.000	12:48:50

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
03/26/2015 12:49:14

50f7
TS4

Standard Lot#: 25414100A4
Last Changed By: J. KNEPPER #46585

WLL

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:49:43
0.100 Cal Check	0.099	12:50:04
Air Blank	0.000	12:50:33

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

60F7
TSH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

WLL
4.20.2015

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3.26.2015 TIME 1229
Intoxilyzer Serial # 80-006048 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
03/26/2015 12:29:39

70F7
B4

Standard Lot#: 25414100A4
Last Changed By: J. KNEPPER #46585

WLL

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:30:08
0.100 Cal Check	0.100	12:30:29
Air Blank	0.000	12:30:59

Preventative Maintenance
Performed -- Timer Reset