



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten initials and date:*  
 11-25-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006048 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-5-2015 Cal 1 0.099 AC Cal 2 0.100 AC

Between the dates of 10-8-2015 and 11-5-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

*Terry Gallegos*

Terry S. Gallegos


11-25-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25<sup>th</sup>, day of November, 2015.

*Elma Escarcega*

Notary Public

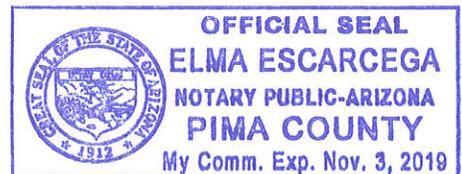

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TS4  
AA  
11-25-15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 11-5-2015 TIME 1342

INTOXILYZER SERIAL # 80-006048 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
11/05/2015

13:45:18

2057  
TSh  
JL

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:46:07
Diagnostic Test	Pass	13:46:42
Air Blank	0.000	13:47:09
0.100 Cal Check	0.100	13:47:30
Air Blank	0.000	13:47:59
Subject Test	0.000	13:48:14
Air Blank	0.000	13:48:42
Wait	ABT*	
Air Blank	0.000	13:49:20

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
11/05/2015

13:49:23

3 of 7  
TS4  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:50:35
Diagnostic Test	Pass	13:51:10
Air Blank	0.000	13:51:38
0.100 Cal Check	RFI*	13:51:55
Air Blank	IPS**	13:52:04

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
11/05/2015

13:52:30

Yof 7  
136  
AA

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:53:17
Diagnostic Test	Pass	13:53:52
Air Blank	0.000	13:54:19
0.100 Cal Check	0.100	13:54:40
Air Blank	0.000	13:55:09
Subject Test	0.079*	13:58:17
Air Blank	0.000	13:58:49
Air Blank	0.000	13:59:16
Subject Test	XXX**	13:59:56
Air Blank	0.000	14:00:30

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
11/05/2015 14:02:46

Soft  
TS4  
JK

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:03:15
0.100 Cal Check	0.100	14:03:36
Air Blank	0.000	14:04:05

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6087  
TS4  
11-25-15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 11-5-2015 TIME 1342  
Intoxilyzer Serial # 80-006048 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C.
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
11/05/2015 13:42:53

70f7  
TS4  


-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

OAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:43:46
0.100 Cal Check	0.099	13:44:07
Air Blank	0.000	13:44:37

Preventative Maintenance  
Performed -- Timer Reset