



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Handwritten initials and date:
 12-18-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006048 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12.7.2015 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 11.5.2015 and 12.7.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Handwritten signature of Terry S. Gallegos

Terry S. Gallegos


12.17.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2015.

Handwritten signature of Elma Escarcega

Notary Public

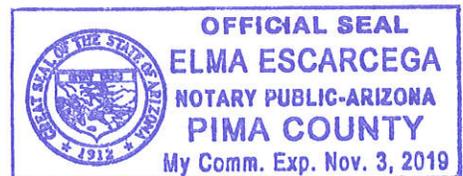

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
TS4
12-18-15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-7-2015 TIME 1046
INTOXILYZER SERIAL # 80-006048 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS > 31 days between checks TS4

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
TSh


INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
12/07/2015

10:52:33

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:53:23 |
| Diagnostic Test | Pass | 10:53:57 |
| Air Blank | 0.000 | 10:54:25 |
| 0.100 Cal Check | 0.102 | 10:54:46 |
| Air Blank | 0.000 | 10:55:15 |
| Subject Test | 0.000 | 10:55:41 |
| Air Blank | 0.000 | 10:56:08 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:56:42 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7
Tg
JL

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
12/07/2015

10:56:48

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:58:11 |
| Diagnostic Test | RFI* | 10:58:22 |
| Air Blank | IPS** | 10:58:34 |

* RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7
TS4
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
12/07/2015

11:00:07

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:00:55 |
| Diagnostic Test | Pass | 11:01:30 |
| Air Blank | 0.000 | 11:01:58 |
| 0.100 Cal Check | 0.101 | 11:02:18 |
| Air Blank | 0.000 | 11:02:48 |
| Subject Test | 0.167* | 11:05:55 |
| Air Blank | 0.000 | 11:06:29 |
| Air Blank | 0.000 | 11:06:56 |
| Subject Test | XXX** | 11:07:27 |
| Air Blank | 0.000 | 11:08:01 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Sof 7
TS4
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
12/07/2015 11:08:25

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:08:55 |
| 0.100 Cal Check | 0.100 | 11:09:16 |
| Air Blank | 0.000 | 11:09:45 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
TSK
12-18-15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-7-2015 TIME 1046
Intoxilyzer Serial # 80 006048 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads ~~"PUSH BUTTON TO START"~~ "Contact QAS 31 Day PREV Maint" TSK
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7
TS4
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006048
Core Version: 8105.48
12/07/2015 10:46:05

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:46:50 |
| 0.100 Cal Check | 0.101 | 10:47:11 |
| Air Blank | 0.000 | 10:47:40 |

Preventative Maintenance
Performed -- Timer Reset