



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Alb*  
1/11/16

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006048 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-30-2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12-7-2015 and 12-30-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —  
The instrument was returned to service on —

*Terry S. Gallegos*  
Terry S. Gallegos


1-7-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 7<sup>th</sup>, day of January, 2016.

*Edma Escarcega*

Notary Public

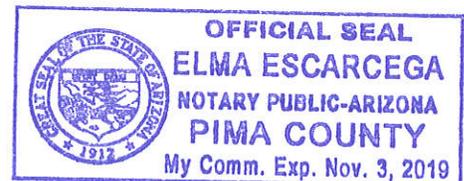

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TSJ  
als  
1/11/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-30-2015 TIME 1245

INTOXILYZER SERIAL # 80-006048 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
12/30/2015

12:47:48

2 of 7  
Tsh  
all

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:48:35
Diagnostic Test	Pass	12:49:09
Air Blank	0.000	12:49:37
0.100 Cal Check	0.100	12:49:58
Air Blank	0.000	12:50:27
Subject Test	0.000	12:51:09
Air Blank	0.000	12:51:36
Wait	ABT*	
Air Blank	0.000	12:52:10

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
12/30/2015

12:53:44

30f7  
Bly  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:54:44
Diagnostic Test	Pass	12:55:19
Air Blank	0.000	12:55:47
0.100 Cal Check	RFI*	12:55:59
Air Blank	IPS**	12:56:14

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
12/30/2015

12:57:09

4 of 7  
By  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:58:11
Diagnostic Test	Pass	12:58:46
Air Blank	0.000	12:59:14
0.100 Cal Check	0.100	12:59:34
Air Blank	0.000	13:00:04
Subject Test	0.038*	13:03:11
Air Blank	0.000	13:03:43
Air Blank	0.000	13:04:10
Subject Test	XXX**	13:04:27
Air Blank	0.000	13:05:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
12/30/2015 13:05:13

*5 of 7*  
*B4*  
*[Signature]*

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:05:42
0.100 Cal Check	0.099	13:06:02
Air Blank	0.000	13:06:32

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60 of 7  
TS4  
AD  
1/11/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 12-30-2015 TIME 1245  
Intoxilyzer Serial # 80-006048 Location DUI-HQ

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (  ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (  ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (  ) 4. Air blank completed.
- (  ) 5. Calibration check completed. Test results 0. 100 AC.
- (  ) 6. Air blank completed.
- (  ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006048  
Core Version: 8105.48  
12/30/2015 12:45:10

7off  
By  
*[Signature]*

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:45:38
0.100 Cal Check	0.100	12:45:59
Air Blank	0.000	12:46:28

Preventative Maintenance  
Performed -- Timer Reset