



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 2.3.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/2/14 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12/5/13 and 1/2/14 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione
 Sarah L. Cione

1/9/14
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th day of January, 2014.

Elma Escarcega
 Notary Public

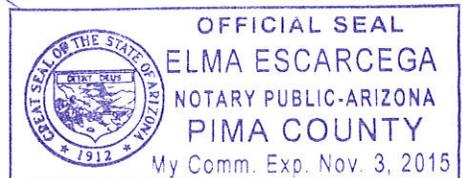

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc
TS4
2.3.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/2/14 TIME 932

INTOXILYZER SERIAL # 80-006050 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
01/02/2014 09:34:35

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:35:18
Diagnostic Test	Pass	09:35:53
Air Blank	0.000	09:36:20
0.100 Cal Check	0.101	09:36:41
Air Blank	0.000	09:37:10
Subject Test	0.000	09:37:27
Air Blank	0.000	09:37:55
Wait	ABT*	
Air Blank	0.000	09:38:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 SC
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
01/02/2014 09:39:59

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:43
Diagnostic Test	Pass	09:41:17
Air Blank	0.000	09:41:45
0.100 Cal Check	RFI*	09:42:00
Air Blank	IPS**	09:42:09

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc
Bk

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
01/02/2014

09:42:15

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:42:57
Diagnostic Test	Pass	09:43:32
Air Blank	0.000	09:43:59
0.100 Cal Check	0.099	09:44:20
Air Blank	0.000	09:44:49
Subject Test	0.630*	09:47:57
Air Blank	0.000	09:48:31
Air Blank	0.000	09:48:58
Subject Test	XXX**	09:49:18
Air Blank	0.000	09:49:52

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
01/02/2014 09:51:01

5 of 7 sc
BSH

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:51:29
0.100 Cal Check	0.099	09:51:50
Air Blank	0.000	09:52:19

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc
Bh
2.3.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/2/14 TIME 932
Intoxilyzer Serial # 80-006050 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
01/02/2014 09:32:54

Bh 7 of 7 sc

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:33:22
0.100 Cal Check	0.100	09:33:43
Air Blank	0.000	09:34:12

Preventative Maintenance
Performed -- Timer Reset