



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS6
 4.14.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/27/14 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2/26/14 and 3/27/14 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione
 Sarah L. Cione

4/1/14
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15th day of April, 2014.

Elma Escarcega
 Notary Public

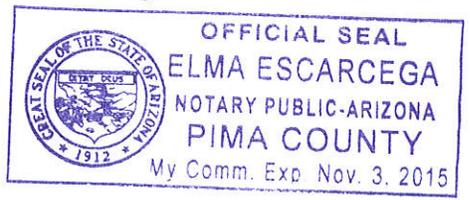

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc
TS4
4.14.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/27/14 TIME 852

INTOXILYZER SERIAL # 80-006050 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
03/27/2014

08:54:33

2 of 7_{sc}
TS4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:55:39
Diagnostic Test	Pass	08:56:14
Air Blank	0.000	08:56:41
0.100 Cal Check	0.100	08:57:02
Air Blank	0.000	08:57:31
Subject Test	0.000	08:58:16
Air Blank	0.000	08:58:43
Wait	ABT*	
Air Blank	0.000	09:00:00

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
03/27/2014

09:00:14

3 of 7 SL
TS4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:01:15
Diagnostic Test	Pass	09:01:50
Air Blank	0.000	09:02:18
0.100 Cal Check	0.100	09:02:38
Air Blank	0.000	09:03:08
Subject Test	RFI*	09:03:33
Air Blank	IPS**	09:03:44

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
03/27/2014

09:14:24

4 of 7_{sc}
TBG

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:15:19
Diagnostic Test	Pass	09:15:54
Air Blank	0.000	09:16:21
0.100 Cal Check	0.100	09:16:42
Air Blank	0.000	09:17:11
Subject Test	0.146*	09:20:19
Air Blank	0.000	09:20:51
Air Blank	0.000	09:21:18
Subject Test	XXX**	09:21:49
Air Blank	0.000	09:22:23

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
03/27/2014 09:22:54

5 of 7 sc
TS4.

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:23:22
0.100 Cal Check	0.099	09:23:43
Air Blank	0.000	09:24:12

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

T34

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

4.14.2014

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/27/14 TIME 852
Intoxilyzer Serial # 80-006050 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
03/27/2014 08:52:06

7 of 7 sc
TSy
4.14.2014

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:52:34
0.100 Cal Check	0.100	08:52:55
Air Blank	0.000	08:53:25

Preventative Maintenance
Performed -- Timer Reset