



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TJG  
 10-29-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/9/14 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 9/12/14 and 10/9/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/13/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of October, 2014.

Elma Escarcega

Notary Public

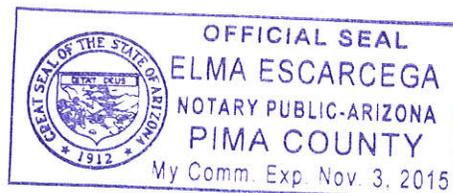

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
sc  
TS4  
10.29.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/9/14 TIME 1035

INTOXILYZER SERIAL # 80-006050 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
10/09/2014

10:37:15

2 of 7  
sc  
154

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:37:57
Diagnostic Test	Pass	10:38:32
Air Blank	0.000	10:39:00
0.100 Cal Check	0.098	10:39:21
Air Blank	0.000	10:39:50
Subject Test	0.000	10:40:25
Air Blank	0.000	10:40:52
Wait	ABT*	
Air Blank	0.000	10:42:07

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
10/09/2014

10:43:34

3 of 7 sc  
sh

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:44:34
Diagnostic Test	Pass	10:45:09
Air Blank	0.000	10:45:36
0.100 Cal Check	0.099	10:45:57
Air Blank	0.000	10:46:27
Subject Test	RFI*	10:46:57
Air Blank	IPS**	10:47:08

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
10/09/2014

10:47:14

4 of 7 sc  
154

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:48:13
Diagnostic Test	Pass	10:48:48
Air Blank	0.000	10:49:16
0.100 Cal Check	0.100	10:49:36
Air Blank	0.000	10:50:06
Subject Test	0.167*	10:53:13
Air Blank	0.000	10:53:46
Air Blank	0.000	10:54:13
Subject Test	XXX**	10:55:36
Air Blank	0.000	10:56:09

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
10/09/2014 10:56:59

5 of 7  
sc  
Bh

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:57:27
0.100 Cal Check	0.099	10:57:48
Air Blank	0.000	10:58:18

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TSh sc  
10-29-2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/9/14 TIME 1035  
Intoxilyzer Serial # 80-006050 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
10/09/2014 10:35:36

7 of 7 SC  
DH

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:36:05
0.100 Cal Check	0.099	10:36:26
Air Blank	0.000	10:36:55

Preventative Maintenance  
Performed -- Timer Reset