



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
 2/16/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-29-2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 1-2-2015 and 1-29-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


2-12-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12<sup>th</sup>, day of February, 20 15.

Elma Escarcega

Notary Public

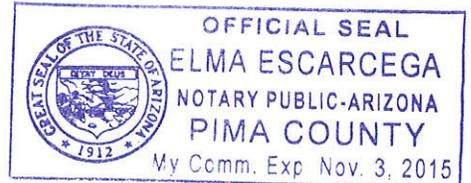

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10F7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
2/16/15

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 1.29.2015 TIME 1420

INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/29/2015

14:23:47

2 of 7  
TS  
SC

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:24:32
Diagnostic Test	Pass	14:25:07
Air Blank	0.000	14:25:35
0.100 Cal Check	0.100	14:25:56
Air Blank	0.000	14:26:25
Subject Test	0.000	14:26:41
Air Blank	0.000	14:27:09
Wait	ABT*	
Air Blank	0.000	14:27:45

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/29/2015

14:29:01

3 of 7  
TJ  
SC

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:29:47
Diagnostic Test	Pass	14:30:22
Air Blank	0.000	14:30:50
0.100 Cal Check	RFI*	14:31:05
Air Blank	IPS**	14:31:13

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7  
Tdy  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/29/2015

14:34:46

---

Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:35:36
Diagnostic Test	Pass	14:36:11
Air Blank	0.000	14:36:38
0.100 Cal Check	0.100	14:36:59
Air Blank	0.000	14:37:29
Subject Test	0.088*	14:40:36
Air Blank	0.000	14:41:08
Air Blank	0.000	14:41:36
Subject Test	XXX**	14:41:55
Air Blank	0.000	14:42:29

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/29/2015 14:43:13

5 of 7  
TJ4  
SC

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:43:42
0.100 Cal Check	0.099	14:44:03
Air Blank	0.000	14:44:32

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
2/16/15

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28843 AGENCY Tucson Police Department  
DATE 1-29-2015 TIME 1420  
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/29/2015 14:20:38

7077  
TS4  
SC

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:21:09
0.100 Cal Check	0.100	14:21:30
Air Blank	0.000	14:22:00

Preventative Maintenance  
Performed -- Timer Reset