



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

Handwritten initials and date: 7/2/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.18.2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 5.22.2015 and 6.18.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Handwritten signature of Terry S. Gallegos

Terry S. Gallegos

6.24.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24th, day of June, 2015.

Handwritten signature of Notary Public

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

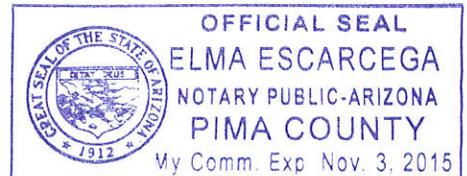


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f7
TSh
Jde
7/2/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 6.18.2015 TIME 1108

INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
06/18/2015

11:10:39

2077
D4
Ad

Standard Lot#: 11514100A2
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:11:26
Diagnostic Test	Pass	11:12:01
Air Blank	0.000	11:12:28
0.100 Cal Check	0.100	11:12:49
Air Blank	0.000	11:13:18
Subject Test	0.000	11:13:34
Air Blank	0.000	11:14:02
Wait	ABT*	
Air Blank	0.000	11:14:35

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
06/18/2015

11:18:22

3 of 7
TJ
AD

Standard Lot#: 11514100A2
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:19:18
Diagnostic Test	Pass	11:19:53
Air Blank	0.000	11:20:20
0.100 Cal Check	RFI*	11:20:35
Air Blank	IPS**	11:20:50

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
06/18/2015

11:21:48

4 of 7
TS4
JL

Standard Lot#: 11514100A2
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:22:36
Diagnostic Test	Pass	11:23:11
Air Blank	0.000	11:23:39
0.100 Cal Check	0.101	11:23:59
Air Blank	0.000	11:24:29
Subject Test	0.097*	11:27:36
Air Blank	0.000	11:28:08
Air Blank	0.000	11:28:36
Subject Test	XXX**	11:28:54
Air Blank	0.000	11:29:27

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
06/18/2015 11:30:09

50f7
Tg
JL

Standard Lot#: 11514100A2
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:30:38
0.100 Cal Check	0.098	11:30:59
Air Blank	0.000	11:31:28

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6007
TJG
7/2/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.18.2015 TIME 1108
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
06/18/2015 11:08:33

Fof7
TSy
JL

Standard Lot#: 11514100A2
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:09:02
0.100 Cal Check	0.100	11:09:22
Air Blank	0.000	11:09:52

Preventative Maintenance
Performed -- Timer Reset