



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

all
 7/30/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-16-2015 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6-18-2015 and 7-16-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Terry Gallegos

Terry S. Gallegos


7-30-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 30th, day of July, 2015.

Elma Escarcega

Notary Public

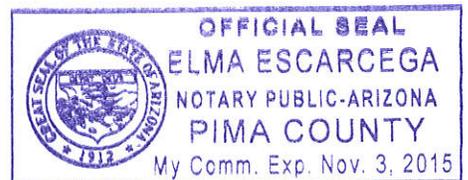

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9
TS
7/13/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 7-16-2015 TIME 1219
INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

20f9
TS4
AD

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/16/2015

12:20:55

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:22:08
Diagnostic Test	Pass	12:22:43
Air Blank	0.000	12:23:10
0.100 Cal Check	0.099	12:23:31
Air Blank	0.000	12:24:01
Subject Test	0.000	12:24:24
Air Blank	0.000	12:24:51
Wait	ABT*	
Air Blank	0.000	12:25:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 9
TS4
JDD

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/16/2015

12:26:01

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:27:47
Diagnostic Test	Pass	12:28:22
Air Blank	0.000	12:28:50
0.100 Cal Check	RFI*	12:29:06
Air Blank	IPS**	12:29:14

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/16/2015

12:29:24

4 of 9
TS4
add

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:30:11
Diagnostic Test	Pass	12:30:46
Air Blank	0.000	12:31:14
0.100 Cal Check	0.099	12:31:35
Air Blank	0.000	12:32:04
Subject Test	0.131*	12:35:11
Air Blank	0.000	12:35:44
Air Blank	0.000	12:36:11
Subject Test	XXX**	12:36:37
Air Blank	0.000	12:37:11

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/16/2015 12:37:41

50f9
TS4
[Signature]

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:38:09
0.100 Cal Check	0.100	12:38:30
Air Blank	0.000	12:39:00

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9
TS4
DLE
7/30/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 7-16-2015 TIME 1219
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- () 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- () 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed.
- () 5. Calibration check completed. Test results 0. 101 AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/16/2015 12:19:17

7059
TSH
Jed

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:19:47
0.100 Cal Check	0.101	12:20:07
Air Blank	0.000	12:20:37

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

80f9
TS4
Alb
7/20/15

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 7/9/15
QAS: HAWKINS 42401
Intoxilyzer Serial Number: 80-006050
New Lot Number: 25414100A4
New Expiration Date: 10/5/16

- (X) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (X) 2. Remove the old gas tank and install the new one.
- (X) 3. Hit the "ESC" key twice and swipe the QAS card.
- (X) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (X) 5. Hit the "ESC" key twice.
- (X) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (X) 7. Air blank completed.
- (X) 8. Calibration check completed. Test results 0. 098 AC.
- (X) 9. Air blank completed.
- (X) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

9 of 9
754
JLB

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
07/09/2015 20:40:39

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: J. HAWKINS #42401
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	20:41:12
0.100 Cal Check	0.098	20:41:33
Air Blank	0.000	20:42:02

JLB
Preventative Maintenance NOT
Performed -- Timer NOT Reset