



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Allo
 12-18-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-3-2015 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 11-5-2015 and 12-3-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


12.17-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2015.

Elma Escarcega

Notary Public

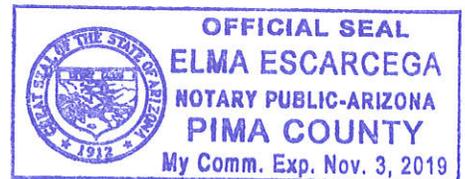

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
TJ
12-18-15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-3-2015 TIME 1211

INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015

12:12:48

*2018
Dh
Ad*

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:13:40
Diagnostic Test	Pass	12:14:18
Air Blank	0.000	12:14:45
0.100 Cal Check	0.100	12:15:06
Air Blank	0.000	12:15:35
Subject Test	0.000	12:16:14
Air Blank	0.000	12:16:41
Wait	ABT*	
Air Blank	0.000	12:17:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015

12:17:20

3058
Tdy
JL

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:18:11
Diagnostic Test	Pass	12:18:46
Air Blank	RFI*	12:18:52
Air Blank	RFI*	12:19:09

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8
D4
ald

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015 13:33:51

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST_SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	13:34:23
Air Blank	0.000	13:34:52

*Improper Sample

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015

12:20:46

5078
T54
Ad

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:21:39
Diagnostic Test	Pass	12:22:14
Air Blank	0.000	12:22:42
0.100 Cal Check	0.102	12:23:02
Air Blank	0.000	12:23:32
Subject Test	0.056*	12:26:39
Air Blank	0.000	12:27:11
Air Blank	0.000	12:27:39
Subject Test	XXX**	12:27:57
Air Blank	0.000	12:28:30

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015 12:28:51

60f8
TJ
[Signature]

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:29:20
0.100 Cal Check	0.101	12:29:41
Air Blank	0.000	12:30:10

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

Zof8
TJ
Ab
12-18-15

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-3-2015 TIME 1211
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/03/2015 12:11:14

80f8
TJ
[Signature]

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:11:43
0.100 Cal Check	0.102	12:12:04
Air Blank	0.000	12:12:33

Preventative Maintenance
Performed -- Timer Reset