



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 1-13-2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/31/15 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 12/3/15 and 12/31/15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

[Handwritten Signature]

Amy L. Autret


1/11/16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 11, day of January, 2016.

[Handwritten Signature]

Notary Public

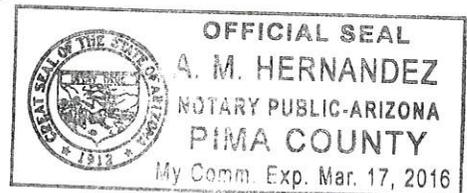

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
JDD
TS4
1-13-2016

QA SPECIALIST Amy Autrot 53546 AGENCY Tucson Police Department

DATE 12-31-15 TIME 1424

INTOXILYZER SERIAL # 80-006050 LOCATION Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

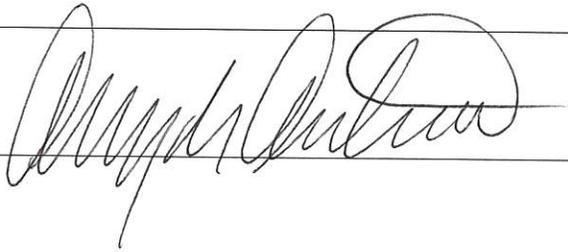

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/31/2015

14:27:20

2 of 7
AS
TS4

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:28:45
Diagnostic Test	Pass	14:29:20
Air Blank	0.000	14:29:47
0.100 Cal Check	0.101	14:30:08
Air Blank	0.000	14:30:38
Subject Test	0.000	14:30:58
Air Blank	0.000	14:31:26
Wait	ABT*	
Air Blank	0.000	14:32:01

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/31/2015

14:32:13

3 of 7
JLH
Bh

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:33:04
Diagnostic Test	Pass	14:33:39
Air Blank	0.000	14:34:06
0.100 Cal Check	RFI*	14:34:20
Air Blank	IPS**	14:34:29

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/31/2015

14:35:17

4 of 7
JL
TS

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:36:21
Diagnostic Test	Pass	14:36:56
Air Blank	0.000	14:37:24
0.100 Cal Check	0.102	14:37:44
Air Blank	0.000	14:38:14
Subject Test	0.103*	14:41:21
Air Blank	0.000	14:41:54
Air Blank	0.000	14:42:21
Subject Test	XXX**	14:42:49
Air Blank	0.000	14:43:23

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/31/2015 14:44:03

5 of 7
JL
TS4

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:44:32
0.100 Cal Check	0.100	14:44:53
Air Blank	0.000	14:45:22

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

leaf 7
TS4
1-13-2016

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 12-31-15 TIME 1424

Intoxilyzer Serial # 80-006050 Location Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist

SIGNATURE



DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006050
Core Version: 8105.48
12/31/2015

14:24:34

Standard Lot#: 25414100A4
Last Changed By: J. HAWKINS #42401

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:25:03
0.100 Cal Check	0.101	14:25:24
Air Blank	0.000	14:25:53

Preventative Maintenance
Performed -- Timer Reset

7 of 7
[Signature]
TS4
1.13.2016