



WLL  
2/11/16

# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-28-2016 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 12-31-2015 and 1-28-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry S. Gallegos

2-9-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 9th, day of February, 2016.

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

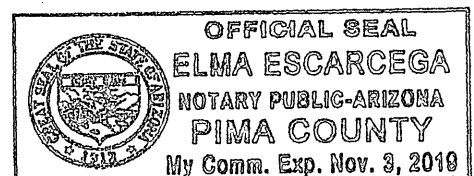


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f8  
TB4  
WLL  
2/11/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 1-28-2016 TIME 1026

INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.  
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.  
(☒) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(☒) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(☒) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016

10:30:43

2 of 8  
Bh  
WLL

Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:31:35
Diagnostic Test	Pass	10:32:10
Air Blank	0.000	10:32:38
0.100 Cal Check	0.101	10:32:58
Air Blank	0.000	10:33:28
Subject Test	0.000	10:34:01
Air Blank	0.000	10:34:28
Wait	ABT*	
Air Blank	0.000	10:35:05

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016

10:46:48

30f8  
BH  
WLL

---

Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:47:41
Diagnostic Test	Pass	10:48:16
Air Blank	0.000	10:48:43
0.100 Cal Check	RFI*	10:48:58
Air Blank	IPS**	10:49:06

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016

10:51:37

40f8  
TS4  
WV

Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:29
Diagnostic Test	Pass	10:53:04
Air Blank	0.000	10:53:32
0.100 Cal Check	0.100	10:53:53
Air Blank	0.000	10:54:22
Subject Test	0.158*	10:57:29
Air Blank	0.000	10:58:01
Air Blank	0.000	10:58:29
Subject Test	DEF**	11:01:36
Air Blank	0.000	11:02:05

\*Deficient Sample

\*\*Deficient Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016

12:37:39

50f8  
TS4  
WLV

---

Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: '111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:39:06
Diagnostic Test	Pass	12:39:41
Air Blank	0.000	12:40:08
0.100 Cal Check	0.100	12:40:29
Air Blank	0.000	12:40:59
Subject Test	XXX*	12:41:25
Air Blank	0.000	12:41:59

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016 12:42:10

60f8  
Bh

WLV

-----  
Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:42:39
0.100 Cal Check	0.099	12:43:00
Air Blank	0.000	12:43:30

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-28-2016 TIME 1026  
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
01/28/2016 10:26:26

-----  
Standard Lot#: 25414100A4  
Last Changed By: J. HAWKINS #42401  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:26:54
0.100 Cal Check	0.101	10:27:15
Air Blank	0.000	10:27:45

Preventative Maintenance  
Performed -- Timer Reset

80f8  
Ph  
WV