



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten initials and date:*  
 9-19-16

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006050 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-8-2016 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 8-11-2016 and 9-8-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

*Handwritten signature of Terry S. Gallegos*

Terry S. Gallegos


9-15-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15<sup>th</sup>, day of September, 2016.

*Handwritten signature of Notary Public*

Notary Public

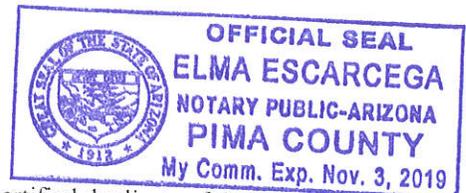

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7  
TS4  
add  
9-19-16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9-8-2016 TIME 1049

INTOXILYZER SERIAL # 80-006050 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 106 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
09/08/2016 10:55:16

20f7  
TS4  
[Signature]

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:07
Diagnostic Test	Pass	10:56:42
Air Blank	0.000	10:57:09
0.100 Cal Check	0.101	10:57:30
Air Blank	0.000	10:57:59
Subject Test	0.000	10:59:06
Air Blank	0.000	10:59:34
Wait	ABT*	
Air Blank	0.000	11:00:08

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
09/08/2016

11:00:30

3 of 7  
TB4  
AD

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:02:13
Diagnostic Test	Pass	11:02:48
Air Blank	0.000	11:03:15
0.100 Cal Check	RFI*	11:03:31
Air Blank	IPS**	11:03:39

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
09/08/2016

11:03:56

4 of 7  
By  
[Signature]

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:04:47
Diagnostic Test	Pass	11:05:21
Air Blank	0.000	11:05:49
0.100 Cal Check	0.101	11:06:10
Air Blank	0.000	11:06:39
Subject Test	0.063*	11:09:46
Air Blank	0.000	11:10:18
Air Blank	0.000	11:10:45
Subject Test	XXX**	11:12:16
Air Blank	0.000	11:12:50

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 6105.46  
09/08/2016 11:45:24

50f7  
734  
*[Signature]*

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:46:04
0.100 Cal Check	0.100	11:46:25
Air Blank	0.000	11:46:54

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
D4  
DLO  
9-19-16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 9.8.2016 TIME 1049  
Intoxilyzer Serial # 80-006050 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

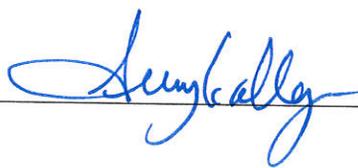SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006050  
Core Version: 8105.48  
09/08/2016 10:49:25

7 of 7  
TS4  
[Signature]

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:49:58
0.100 Cal Check	0.101	10:50:19
Air Blank	0.000	10:50:49

Preventative Maintenance  
Performed -- Timer Reset