



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
12.3.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/7/13 Cal 1 0.100 AC Cal 2 0.101 AC

Between the dates of 10/9/13 and 11/7/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


11/20/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20<sup>th</sup>, day of November, 2013.

Elma Escarcega

Notary Public

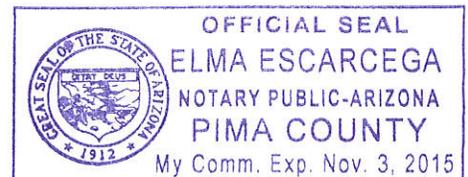

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
sc  
TSy  
12-3-2013

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/7/13 TIME 1105

INTOXILYZER SERIAL # 80-006051 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
11/07/2013

11:11:49

2 of 7  
TS4 SC

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:12:40
Diagnostic Test	Pass	11:13:15
Air Blank	0.000	11:13:42
0.100 Cal Check	0.102	11:14:03
Air Blank	0.000	11:14:32
Subject Test	0.000	11:14:58
Air Blank	0.000	11:15:26
Wait	ABT*	
Air Blank	0.000	11:16:18

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
11/07/2013

11:16:32

3 of 7  
TS4 SC

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:17:25
Diagnostic Test	Pass	11:18:00
Air Blank	0.000	11:18:28
0.100 Cal Check	0.100	11:18:48
Air Blank	0.000	11:19:18
Subject Test	RFI*	11:19:39
Air Blank	IPS**	11:19:47

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-00605i  
Core Version: 8105.48  
11/07/2013

11:25:08

4 of 7  
TS4 SC

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:26:01
Diagnostic Test	Pass	11:26:36
Air Blank	0.000	11:27:03
0.100 Cal Check	0.101	11:27:24
Air Blank	0.000	11:27:54
Subject Test	0.105*	11:31:01
Air Blank	0.000	11:31:33
Air Blank	0.000	11:32:00
Subject Test	XXX**	11:32:25
Air Blank	0.000	11:32:58

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
11/07/2013 11:33:59

5 of 7  
sc  
TB4

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:34:27
0.100 Cal Check	0.101	11:34:48
Air Blank	0.000	11:35:17

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7<sub>SC</sub>  
TBH  
12.3.2013

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/7/13 TIME 1105  
Intoxilyzer Serial # 80-006051 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
11/07/2013 11:05:35

7 of 7 sc  
TSh

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:06:03
0.100 Cal Check	0.100	11:06:24
Air Blank	0.000	11:06:53

Preventative Maintenance  
Performed -- Timer Reset