



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TG  
2.3.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/2/14 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 12/5/13 and 1/2/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/9/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9<sup>th</sup> day of January, 20 14.

Elma Escarcega

Notary Public

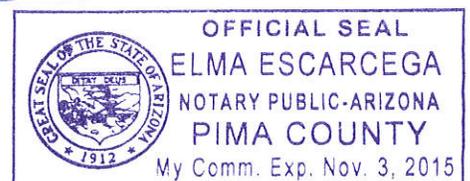

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 sc  
B4  
2.3.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/2/14 TIME 841

INTOXILYZER SERIAL # 80-006051 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(  ) 1. Clock time check.

(  ) 2. Date check.

OPERATIONAL TESTS

(  ) 1. Alcohol-free subject test result 0. 000 AC.

(  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
01/02/2014

08:44:17

2 of 7<sub>sc</sub>  
134

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:45:00
Diagnostic Test	Pass	08:45:35
Air Blank	0.000	08:46:02
0.100 Cal Check	0.100	08:46:23
Air Blank	0.000	08:46:52
Subject Test	0.000	08:47:35
Air Blank	0.000	08:48:03
Wait	ABT*	
Air Blank	0.000	08:48:38

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
01/02/2014

08:51:49

3 of 7 sc  
Bh

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:52:51
Diagnostic Test	Pass	08:53:26
Air Blank	0.000	08:53:53
0.100 Cal Check	0.100	08:54:14
Air Blank	0.000	08:54:43
Subject Test	RFI*	08:55:00
Air Blank	IPS**	08:55:07

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
01/02/2014

08:56:05

4 of 7 SC  
B4

---

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:57:04
Diagnostic Test	Pass	08:57:39
Air Blank	0.000	08:58:06
0.100 Cal Check	0.101	08:58:27
Air Blank	0.000	08:58:57
Subject Test	0.219*	09:02:04
Air Blank	0.000	09:02:36
Air Blank	0.000	09:03:04
Subject Test	XXX**	09:03:28
Air Blank	0.000	09:04:03

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
01/02/2014 09:04:33

5 of 7 sc  
T84

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:05:01
0.100 Cal Check	0.099	09:05:22
Air Blank	0.000	09:05:51

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
2.3.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/2/14 TIME 841  
Intoxilyzer Serial # 80-006051 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
01/02/2014 08:41:14

7 of 7 SC  
15h

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:41:44
0.100 Cal Check	0.101	08:42:05
Air Blank	0.000	08:42:35

Preventative Maintenance  
Performed -- Timer Reset