



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
12/19/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12.4.2014 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 11.6.2014 and 12.4.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


12-17-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of December, 2014.

Elma Escarcega

Notary Public

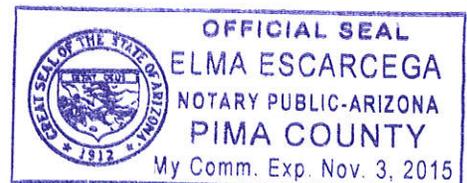

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10f9  
B4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
12/19/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-4-2014 TIME 1554

INTOXILYZER SERIAL # 80-006051 LOCATION ODW-DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE 

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014

15:56:12

2 of 9  
TD  
5c

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 15:57:53 |
| Diagnostic Test | Pass   | 15:58:27 |
| Air Blank       | 0.000  | 15:58:55 |
| 0.100 Cal Check | 0.099  | 15:59:16 |
| Air Blank       | 0.000  | 15:59:45 |
| Subject Test    | 0.000  | 16:00:09 |
| Air Blank       | 0.000  | 16:00:36 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 16:01:10 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014

16:01:17

30 of 9  
TS4

5c

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 16:03:09 |
| Diagnostic Test | Pass   | 16:03:44 |
| Air Blank       | 0.000  | 16:04:12 |
| 0.100 Cal Check | RFI*   | 16:04:27 |
| Air Blank       | IPS**  | 16:04:35 |

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014

16:04:54

4 of 9  
D4

5c

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 16:05:46 |
| Diagnostic Test | Pass   | 16:06:21 |
| Air Blank       | 0.000  | 16:06:49 |
| 0.100 Cal Check | 0.099  | 16:07:09 |
| Air Blank       | 0.000  | 16:07:39 |
| Subject Test    | 0.156* | 16:10:46 |
| Air Blank       | 0.000  | 16:11:20 |
| Air Blank       | 0.000  | 16:11:48 |
| Subject Test    | XXX**  | 16:12:16 |
| Air Blank       | 0.000  | 16:12:51 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014 16:13:12

50f9  
Dy

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 16:13:42 |
| 0.100 Cal Check | 0.098  | 16:14:02 |
| Air Blank       | 0.000  | 16:14:32 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

60f9  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
12/19/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 12.4.2014 TIME 1554  
Intoxilyzer Serial # 80-006051 Location DDW-DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014 15:54:14

70f9  
D4  
SC

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 15:54:44 |
| 0.100 Cal Check | 0.098  | 15:55:04 |
| Air Blank       | 0.000  | 15:55:34 |

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8059  
TD4

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

SC  
12/19/14

Date: 12.4.2014  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-006051  
New Lot Number: 11514100A2  
New Expiration Date: 05.01.2016

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006051  
Core Version: 8105.48  
12/04/2014 16:46:23

9 of 9  
TSh

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

Sc

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 16:46:53 |
| 0.100 Cal Check | 0.100  | 16:47:14 |
| Air Blank       | 0.000  | 16:47:43 |

Preventative Maintenance  
Performed -- Timer Reset