



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

SC
 2/16/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1.29.2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12.31.2014 and 1.29.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos


2.12.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12th, day of February, 20 15.

Elma Escarcega

Notary Public

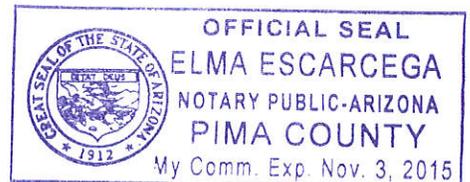

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

10 of 7
TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 3000

SC
2/16/15

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-29-2015 TIME 1517
INTOXILYZER SERIAL # 80-006051 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
01/29/2015

15:19:49

2 of 7
Bg
5c

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:20:40
Diagnostic Test	Pass	15:21:15
Air Blank	0.000	15:21:42
0.100 Cal Check	0.100	15:22:03
Air Blank	0.000	15:22:33
Subject Test	0.000	15:23:11
Air Blank	0.000	15:23:39
Wait	ABT*	
Air Blank	0.000	15:24:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
01/29/2015

15:25:31

30f7
TS4
5C

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:27:15
Diagnostic Test	Pass	15:27:50
Air Blank	0.000	15:28:17
0.100 Cal Check	RFI*	15:28:34
Air Blank	IPS**	15:28:48

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
01/29/2015

15:29:34

4 of 7
D4
SC

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:30:23
Diagnostic Test	Pass	15:30:58
Air Blank	0.000	15:31:25
0.100 Cal Check	0.101	15:31:46
Air Blank	0.000	15:32:15
Subject Test	0.346*	15:35:23
Air Blank	0.000	15:36:04
Air Blank	0.000	15:36:32
Subject Test	XXX**	15:36:59
Air Blank	0.000	15:37:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
01/29/2015 15:38:44

5 of 7
B4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:39:12
0.100 Cal Check	0.099	15:39:33
Air Blank	0.000	15:40:03

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000


STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TS4

SC
2/16/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-29-2015 TIME 1517
Intoxilyzer Serial # 80-006051 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

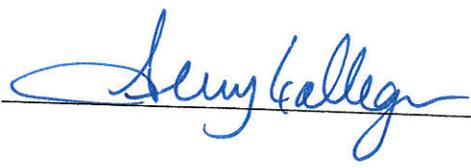SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
01/29/2015 15:17:20

70f7
T6

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:17:49
0.100 Cal Check	0.100	15:18:10
Air Blank	0.000	15:18:40

Preventative Maintenance
Performed -- Timer Reset