



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

AB
 10/23/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10.8.2015 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 9.10.2015 and 10.8.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
 The instrument was returned to service on .

Terry Gallegos

Terry S. Gallegos


10.16.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16th, day of October, 2015.

Elma Escarcega

Notary Public

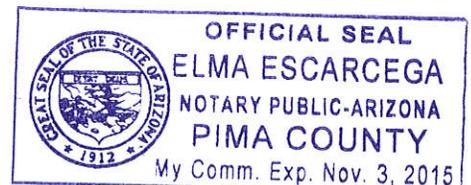

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS
AS
10/23/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10-8-2015 TIME 1514

INTOXILYZER SERIAL # 80-006051 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
10/08/2015

15:20:30

20 of 7
TS
dd

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:21:15
Diagnostic Test	Pass	15:21:50
Air Blank	0.000	15:22:17
0.100 Cal Check	0.099	15:22:38
Air Blank	0.000	15:23:08
Subject Test	0.000	15:23:39
Air Blank	0.000	15:24:06
Wait	ABT*	
Air Blank	0.000	15:24:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
10/08/2015

15:25:30

30F7
TS4
Ab

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:26:50
Diagnostic Test	Pass	15:27:25
Air Blank	0.000	15:27:53
0.100 Cal Check	RFI*	15:28:07
Air Blank	IPS**	15:28:21

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
10/08/2015

15:29:01

4 of 7
TS4
AD

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:29:46
Diagnostic Test	Pass	15:30:21
Air Blank	0.000	15:30:49
0.100 Cal Check	0.098	15:31:09
Air Blank,	0.000	15:31:39
Subject Test	0.081*	15:34:46
Air Blank	0.000	15:35:20
Air Blank	0.000	15:35:47
Subject Test	XXX**	15:36:06
Air Blank	0.000	15:36:39

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
10/08/2015 15:36:52

50f7
TS4
dd

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:37:21
0.100 Cal Check	0.098	15:37:42
Air Blank	0.000	15:38:11

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60f7
TS4
All
10/23/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 10-8-2015 TIME 1514
Intoxilyzer Serial # 80-006051 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
10/08/2015 15:14:00

7 of 7
TS4
JD

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:14:30
0.100 Cal Check	0.098	15:14:51
Air Blank	0.000	15:15:21

Preventative Maintenance
Performed -- Timer Reset