



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

NO
 11-25-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006051 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-5-2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 10-8-2015 and 11-5-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.


Terry S. Gallegos

Terry S. Gallegos

11-25-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25th day of November, 20 15.

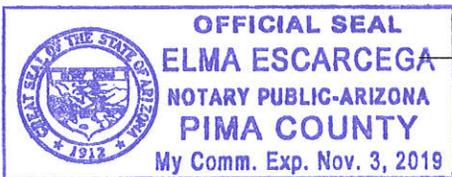

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10 of 10
TSg
[Signature]

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 11-5-2015 TIME 1451
INTOXILYZER SERIAL # 80-006051 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS 10-15-2015:
Officer Szelewski sent a message that the
cal check measured 0.113 AC. Contact QAS message
on Instrument. TSg

SIGNATURE [Signature]

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

20810
TJ
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 15:55:08

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:56:02
Diagnostic Test	Pass	15:56:39
Air Blank	0.000	15:57:06
0.100 Cal Check	0.099	15:57:27
Air Blank	0.000	15:57:57
Subject Test	0.000	15:58:19
Air Blank	0.000	15:58:47
Wait	ABT*	
Air Blank	0.000	15:59:20

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 10
T34
JL

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 15:59:49

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:02:54
Diagnostic Test	Pass	16:03:32
Air Blank	0.000	16:03:59
0.100 Cal Check	RFI*	16:04:13
Air Blank	IPS**	16:04:22

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 10
TS4
DWB

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 16:11:41

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:12:28
Diagnostic Test	Pass	16:13:06
Air Blank	0.000	16:13:33
0.100 Cal Check	0.100	16:13:54
Air Blank	0.000	16:14:23
Subject Test	0.049*	16:17:31
Air Blank	0.000	16:18:04
Air Blank	0.000	16:18:31
Subject Test	0.239**	16:18:59
Air Blank	0.000	16:19:33
Wait		
Air Blank	0.000	16:24:27
Subject Test	0.307**	16:24:57
Air Blank	0.000	16:25:35

*Deficient Sample
**No .020 agreement

Not a Successfully
Completed
Test Sequence

*I ran out of breath
On Invalid sample attempt. Tried again with same
result - I ran out of breath. TS4*

5 of 10
TS4
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 16:31:13

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST_SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:32:05
Diagnostic Test	Pass	16:32:42
Air Blank	0.000	16:33:09
0.100 Cal Check	0.100	16:33:30
Air Blank	0.000	16:34:00
Subject Test	0.000	16:34:41
Air Blank	0.000	16:35:08
Wait		
Air Blank	0.000	16:40:09
Subject Test	XXX*	16:40:31
Air Blank	0.000	16:41:07

*Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

*Subject Test = 0.000
Wait 5 min after Invalid sample attempt,
then did another Invalid Sample TS4.*

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 16:41:43

6 of 10
TJ
JK

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:42:12
0.100 Cal Check	0.100	16:42:33
Air Blank	0.000	16:43:03

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

*Tof 10
TS4
Old
11-25-15*

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 11-5-2015 TIME 1458
Intoxilyzer Serial # 80-006051 Location DUT-HQ

- () 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- () 2. Intoxilyzer 8000 display reads "Contact QAS Cal Check Tolerance" ~~PUSH BUTTON TO START~~ *TS4*
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed. *I used QAS Card to update Std. Then went to Ready for next test. TS4*
- () 5. Calibration check completed. Test results 0. 099 AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8 of 10
DJ
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 14:58:47

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

DAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:00:20
0.100 Cal Check	0.099	15:00:40
Air Blank	0.000	15:01:10

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

9 of 10
TJ
11-25-15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 11-5-2015 TIME 1451
Intoxilyzer Serial # 80006051 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- (✓) 2. Intoxilyzer 8000 display reads "~~PUSH BUTTON TO START~~" ^{"Contact QAS Cal/Check Tolerance"}
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER. ^{TJG}
- (✓) 4. Air blank completed. ^{Used Technician Card, It Did not update std. TJG}
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Instrument Initialization
13:58:28 11/05/2015

10 of 10
TS4
[Signature]

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006051
Core Version: 8105.48
11/05/2015 14:51:42

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	14:52:09
0.100 Cal Check	0.098	14:52:30
Air Blank	0.000	14:53:00