



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
 12/10/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11.13.2014 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 10.15.2014 and 11.13.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.


Terry S. Gallegos

Terry S. Gallegos

12.4.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24th, day of December, 2014.

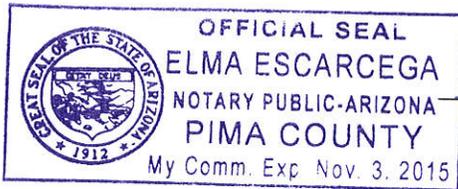

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1057  
TS4  
SC  
12/10/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 2884 AGENCY Tucson Police Department

DATE 11-13-2014 TIME 1028

INTOXILYZER SERIAL # 80-006052 LOCATION ODE - Intox Room

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

2 of 7  
TS4

SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
11/13/2014

10:30:28

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:31:21 |
| Diagnostic Test | Pass   | 10:31:58 |
| Air Blank       | 0.000  | 10:32:25 |
| 0.100 Cal Check | 0.099  | 10:32:46 |
| Air Blank       | 0.000  | 10:33:16 |
| Subject Test    | 0.000  | 10:33:40 |
| Air Blank       | 0.000  | 10:34:07 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 10:34:42 |

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

30F7  
T24

5C

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
11/13/2014 10:36:32

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

OAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:37:24 |
| Diagnostic Test | Pass   | 10:38:01 |
| Air Blank       | 0.000  | 10:38:28 |
| 0.100 Cal Check | 0.099  | 10:38:49 |
| Air Blank       | RFI*   | 10:38:55 |
| Air Blank       | IPS**  | 10:39:03 |

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

40F7  
TSH

SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
11/13/2014 10:39:14

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:40:05 |
| Diagnostic Test | Pass   | 10:40:43 |
| Air Blank       | 0.000  | 10:41:10 |
| 0.100 Cal Check | 0.098  | 10:41:31 |
| Air Blank       | 0.000  | 10:42:00 |
| Subject Test    | 0.012* | 10:45:08 |
| Air Blank       | 0.000  | 10:45:42 |
| Air Blank       | 0.000  | 10:46:09 |
| Subject Test    | XXX**  | 10:46:46 |
| Air Blank       | 0.000  | 10:47:20 |

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
11/13/2014 10:47:43

5C

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:48:12 |
| 0.100 Cal Check | 0.099  | 10:48:32 |
| Air Blank       | 0.000  | 10:49:02 |

Preventative Maintenance  
Performed -- Timer Reset

|

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

60 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

5C  
12/10/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 11-13-2014 TIME 1028  
Intoxilyzer Serial # 80 006052 Location ODE - Intox Room

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

7077  
736

50

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
11/13/2014 10:28:47

Standard Lot #: 25913100A2  
Last Changed By: S. CIONE #50616

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:29:17 |
| 0.100 Cal Check | 0.099  | 10:29:38 |
| Air Blank       | 0.000  | 10:30:07 |

Preventative Maintenance  
Performed -- Timer Reset