



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TJH
1.21.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/9/15 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 12/11/14 and 1/9/15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

[Signature]
Amy L. Autret


1/12/15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 12th, day of January, 2015.

[Signature]

Notary Public

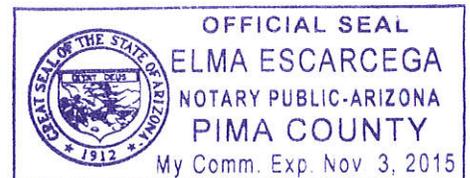

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7
llb
Bh
1-21-2015

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 3000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 1-9-15 TIME 1029

INTOXILYZER SERIAL # 80-006052 LOCATION ODE - Intox Room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

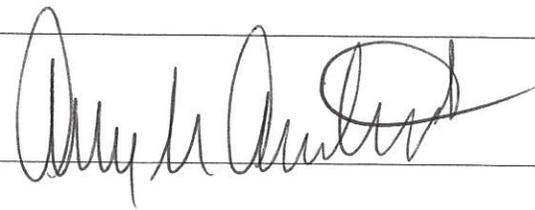SIGNATURE 

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

2 of 7
JL
Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
01/09/2015 10:31:31

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1971
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:32:33
Diagnostic Test	Pass	10:33:07
Air Blank	0.000	10:33:35
0.100 Cal Check	0.100	10:33:56
Air Blank	0.000	10:34:25
Subject Test	0.000	10:34:41
Air Blank	0.000	10:35:08
Wait	ABT*	
Air Blank	0.000	10:35:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7
dlb
Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
01/09/2015 10:36:11

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1971
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:37:14
Diagnostic Test	Pass	10:37:49
Air Blank	0.000	10:38:16
0.100 Cal Check	RFI*	10:38:32
Air Blank	IPS**	10:38:42

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7
JL
TS

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
01/09/2015 10:39:12

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1971
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:40:13
Diagnostic Test	Pass	10:40:48
Air Blank	0.000	10:41:15
0.100 Cal Check	0.099	10:41:36
Air Blank	0.000	10:42:06
Subject Test	0.110*	10:45:13
Air Blank	0.000	10:45:45
Air Blank	0.000	10:46:12
Subject Test	XXX**	10:46:29
Air Blank	0.000	10:47:03

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7
All
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
01/09/2015 10:48:04

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:48:35
0.100 Cal Check	0.099	10:48:56
Air Blank	0.000	10:49:26

Preventative Maintenance
Performed -- Timer Reset

Handwritten initials

*TSK
1-21-2015*

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Ruetz 53546 AGENCY Tucson Police Department

DATE 1-9-15 TIME 1029

Intoxilyzer Serial # 80-006052 Location ODE - Intox Room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Handwritten signature of Amy Ruetz

7 of 7

Handwritten signature

Handwritten initials

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
01/09/2015 10:29:27

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:29:58
0.100 Cal Check	0.099	10:30:19
Air Blank	0.000	10:30:48

Preventative Maintenance
Performed -- Timer Reset