



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 2.11.2015


I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-6-2015 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 1-9-2015 and 2-6-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

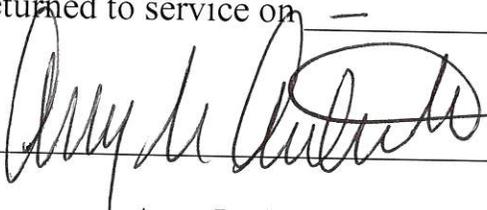

Amy L. Autret


2-10-2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup>, day of February, 20 15.



Notary Public

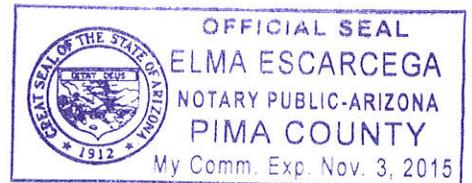

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7  
JLH  
TSY  
2.11.2015

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 2-6-2015 TIME 1156

INTOXILYZER SERIAL # 80-006052 LOCATION ODE - Intox room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS Attempted Invalid Sample twice.

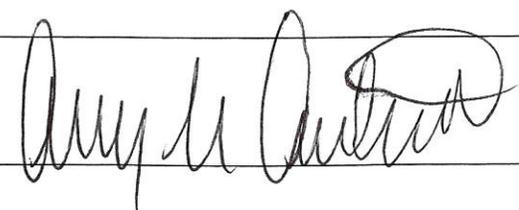SIGNATURE 

2 of 7  
dub  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
02/06/2015 11:59:14

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CLONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/12/1987  
Sex: F Weight: 135  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:00:15
Diagnostic Test	Pass	12:00:52
Air Blank	0.000	12:01:19
0.100 Cal Check	0.100	12:01:40
Air Blank	0.000	12:02:10
Subject Test	0.000	12:02:29
Air Blank	0.000	12:02:57
Wait	ABT*	
Air Blank	0.000	12:03:32

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7  
old  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
02/06/2015 12:04:14

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

OAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/12/1987  
Sex: F Weight: 135  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:05:14
Diagnostic Test	Pass	12:05:51
Air Blank	0.000	12:06:19
0.100 Cal Check	RFI*	12:06:32
Air Blank	IPS**	12:06:42

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 7  
add  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
02/06/2015 12:07:27

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

DAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/12/1987  
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:08:24
Diagnostic Test	Pass	12:09:01
Air Blank	0.000	12:09:29
0.100 Cal Check	0.099	12:09:49
Air Blank	0.000	12:10:19
Subject Test	0.060*	12:13:26
Air Blank	0.000	12:13:58
Air Blank	0.000	12:14:26
Subject Test	0.056	12:15:02
Air Blank	0.000	12:15:34
Wait		
Air Blank	0.000	12:20:30
Subject Test	XXX**	12:20:54
Air Blank	0.000	12:21:27

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

Ran out of breath on  
Invalid sample attempt.  
Waited 5 minutes and  
used breath spray for  
2nd Invalid sample  
attempt.

5 of 7  
all  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
02/06/2015 12:22:41

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:23:13
0.100 Cal Check	0.100	12:23:34
Air Blank	0.000	12:24:03

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE


6 of 7  
JTB  
2.11.2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 2-6-2015 TIME 1156

Intoxilyzer Serial # 80-006052 Location ODE-Intox room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

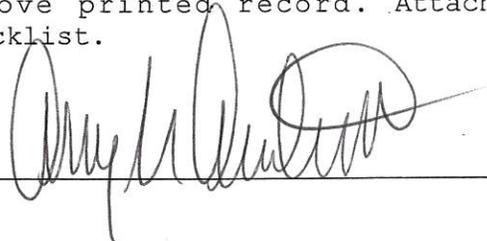SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

7027  
JLB  
734

INTOXILYZER 8000  
Instrument Initialization  
11:46:36 02/06/2015

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
02/06/2015 11:56:42

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:57:14
0.100 Cal Check	0.100	11:57:34
Air Blank	0.000	11:58:04

Preventative Maintenance  
Performed -- Timer Reset