



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
5.7-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-1-2015 Cal 1 0.099 AC Cal 2 0.100 AC

Between the dates of 4-2-2015 and 5-1-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .


Amy L. Autret

5-5-2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 5th, day of May, 2015.

Notary Public

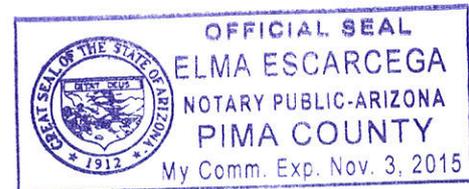

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10F7
dls
TS4
5.7.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 5-1-2015 TIME 0956

INTOXILYZER SERIAL # 80-006052 LOCATION ODE Intox Room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

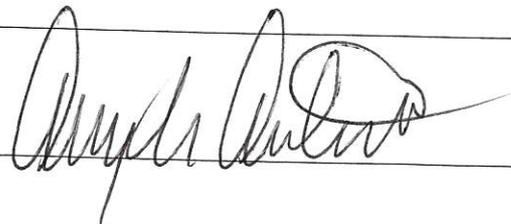

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
MA
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
05/01/2015 09:59:40

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

OAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:00:35
Diagnostic Test	Pass	10:01:13
Air Blank	0.000	10:01:40
0.100 Cal Check	0.100	10:02:01
Air Blank	0.000	10:02:30
Subject Test	0.000	10:02:51
Air Blank	0.000	10:03:18
Wait	ABT*	
Air Blank	0.000	10:03:54

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7
ML

TSy

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
05/01/2015 10:04:35

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:05:33
Diagnostic Test	Pass	10:06:10
Air Blank	0.000	10:06:37
0.100 Cal Check	RFI*	10:06:52
Air Blank	IPS**	10:07:01

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7
Ad
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
05/01/2015 10:07:38

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:08:36
Diagnostic Test	Pass	10:09:13
Air Blank	0.000	10:09:40
0.100 Cal Check	0.101	10:10:01
Air Blank	0.000	10:10:31
Subject Test	0.155*	10:13:38
Air Blank	0.000	10:14:11
Air Blank	0.000	10:14:39
Subject Test	XXX**	10:15:03
Air Blank	0.000	10:15:37

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7
MS
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
05/01/2015 10:16:30

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:17:00
0.100 Cal Check	0.100	10:17:21
Air Blank	0.000	10:17:50

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
MB
By
5.7.2015


QA Specialist Amy Rutret 53546 AGENCY Tucson Police Department

DATE 5-1-2015 TIME 0956

Intoxilyzer Serial # 80-006052 Location ODE Intox Room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

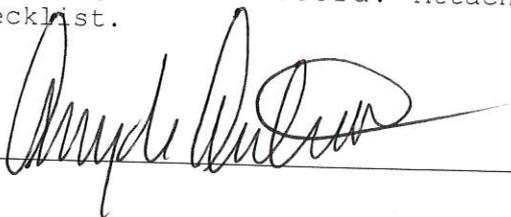

DPS Form Exh G-2 (Rev. 05-01)

7067
JL
TS4

INTOXILYZER 8000
Instrument Initialization
09:13:53 05/01/2015

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
05/01/2015 09:56:26

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

DAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:57:00
0.100 Cal Check	0.099	09:57:21
Air Blank	0.000	09:57:51

Preventative Maintenance
Performed -- Timer Reset