



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
11-3-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/15/2015 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 9/17/2015 and 10/15/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


10/22/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 22nd, day of October, 2015.

Elma Escarcega

Notary Public

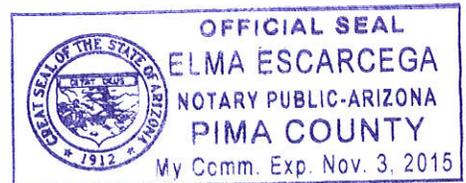

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9
TS4
11-3-2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 10/15/2015 TIME 09:17

INTOXILYZER SERIAL # 80-006052 LOCATION ODE-SGT'S OFFICE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 060 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
10/15/2015 09:20:42

2 of 9
mlv
TS4

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:22:00
Diagnostic Test	Pass	09:22:38
Air Blank	0.000	09:23:05
0.100 Cal Check	0.102	09:23:26
Air Blank	0.000	09:23:56
Subject Test	0.000	09:24:22
Air Blank	0.000	09:24:49
Wait	ABT*	
Air Blank	0.000	09:25:34

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 9
wll
154

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
10/15/2015 09:28:43

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:30:10
Diagnostic Test	Pass	09:30:48
Air Blank	0.000	09:31:15
0.100 Cal Check	RFI*	09:31:30
Air Blank	IPS**	09:31:40

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 9
will
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
10/15/2015 09:33:07

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:34:30
Diagnostic Test	Pass	09:35:08
Air Blank	0.000	09:35:35
0.100 Cal Check	0.100	09:35:56
Air Blank	0.000	09:36:25
Subject Test	0.129*	09:39:33
Air Blank	0.000	09:40:06
Air Blank	0.000	09:40:34
Subject Test	XXX**	09:40:56
Air Blank	0.000	09:41:29

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 9
WU
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
10/15/2015 09:43:30

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:44:00
0.100 Cal Check	0.100	09:44:20
Air Blank	0.000	09:44:50

Preventative Maintenance
Performed -- Timer Reset

|

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 9
WLL
TSy
11.3.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 10/15/2015 TIME 09:17
Intoxilyzer Serial # 80-006052 Location ODE - SGT'S OFFICE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 9
Owll

TS4
11-3-2015

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
10/15/2015 09:17:56

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:18:26
0.100 Cal Check	0.101	09:18:47
Air Blank	0.000	09:19:17

Preventative Maintenance
Performed -- Timer Reset

|



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2
JLB
8 of 9
0 WLL

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TSG
11-3-2015

Date: 9-28-2015

QAS: Amy Autret 53546

Intoxilyzer Serial Number: 80-006052

New Lot Number: 12415100A4

New Expiration Date: 6/5/2017

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:

Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."

- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

2 of 2
AB
9 of 9
WU
TJ

INTOXILYZER 8000
Instrument Initialization
08:29:27 09/28/2015

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006052
Core Version: 8105.48
09/28/2015 08:59:16

Standard Lot#: 12415100A4
Last Changed By: A. AUTRET #53546

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:59:46
0.100 Cal Check	0.100	09:00:07
Air Blank	0.000	09:00:37

Preventative Maintenance NOT
Performed -- Timer NOT Reset