



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 1.13.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006052 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/8/16 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of 12/14/15 and 1/8/16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
 The instrument was returned to service on                     .

*[Handwritten Signature]*

Amy L. Autret


1/11/16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 11, day of January, 2016.

*[Handwritten Signature]*

Notary Public

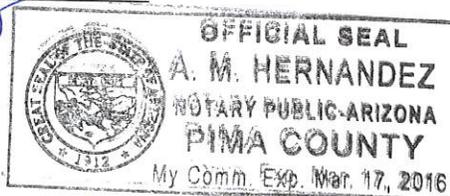

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1047  
JLB  
TS4  
1-13-2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 1/8/16 TIME 0849

INTOXILYZER SERIAL # 80-006052 LOCATION ODE - Intox room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE



Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
01/08/2016

08:51:11

2 of 7  
JWA  
134

Standard Lot#: 12415100A4  
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:52:03
Diagnostic Test	Pass	08:52:38
Air Blank	0.000	08:53:05
0.100 Cal Check	0.102	08:53:26
Air Blank	0.000	08:53:56
Subject Test	0.000	08:54:12
Air Blank	0.000	08:54:39
Wait	ABT*	
Air Blank	0.000	08:55:14

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
01/08/2016

08:55:19

3 of 7  
JLS  
T84

Standard Lot#: 12415100A4  
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:56:13
Diagnostic Test	Pass	08:56:48
Air Blank	0.000	08:57:15
0.100 Cal Check	RFI*	08:57:29
Air Blank	IPS**	08:57:37

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
01/08/2016

08:57:55

4 of 7  
old  
TSY

Standard Lot#: 12415100A4  
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:58:54
Diagnostic Test	Pass	08:59:29
Air Blank	0.000	08:59:56
0.100 Cal Check	0.101	09:00:17
Air Blank	0.000	09:00:47
Subject Test	0.214*	09:03:54
Air Blank	0.000	09:04:27
Air Blank	0.000	09:04:55
Subject Test	XXX**	09:05:16
Air Blank	0.000	09:05:50

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
01/08/2016 09:06:00

507  
JL  
TJ

-----  
Standard Lot#: 12415100A4  
Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:06:30
0.100 Cal Check	0.101	09:06:51
Air Blank	0.000	09:07:21

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE


6 of 7  
JLS  
TS4  
1.13.2016

QA Specialist Amy Antret 53546 AGENCY Tucson Police Department

DATE 1/8/16 TIME 0849

Intoxilyzer Serial # 80-006052 Location ODE - Intox room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

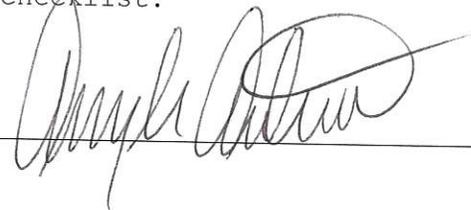SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006052  
Core Version: 8105.48  
01/08/2016 08:49:34

7 of 7  
JLB  
TS4

-----  
Standard Lot#: 12415100A4  
Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:50:04
0.100 Cal Check	0.103	08:50:25
Air Blank	0.000	08:50:54

Preventative Maintenance  
Performed -- Timer Reset