



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

By
2.12.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006053 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/30/14 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 1/2/14 and 1/30/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/11/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11th, day of February, 2014.

Elma Escarcega

Notary Public

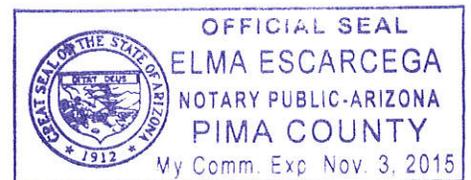

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7_{sc}
TBH
2.12.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/30/14 TIME 837

INTOXILYZER SERIAL # 80-006053 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006053
Core Version: 8105.48
01/30/2014

08:40:41

2 of 7 sc
TB4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:41:40 |
| Diagnostic Test | Pass | 08:42:15 |
| Air Blank | 0.000 | 08:42:42 |
| 0.100 Cal Check | 0.099 | 08:43:03 |
| Air Blank | 0.000 | 08:43:32 |
| Subject Test | 0.000 | 08:43:49 |
| Air Blank | 0.000 | 08:44:17 |
| Wait | ABT* | |
| Air Blank | 0.000 | 08:44:51 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006053
Core Version: 8105.48
01/30/2014

08:48:56

3 of 7 sc
134

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:49:48 |
| Diagnostic Test | Pass | 08:50:23 |
| Air Blank | 0.000 | 08:50:50 |
| 0.100 Cal Check | 0.099 | 08:51:11 |
| Air Blank | 0.000 | 08:51:41 |
| Subject Test | RFI* | 08:52:32 |
| Air Blank | IPS** | 08:52:42 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006053
Core Version: 8105.48
01/30/2014

08:53:32

4 of 7 sc
TSK

Standard Lot#: 19613100A4 #50616
Last Changed By: S. CIONE

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:54:17 |
| Diagnostic Test | Pass | 08:54:52 |
| Air Blank | 0.000 | 08:55:19 |
| 0.100 Cal Check | 0.100 | 08:55:40 |
| Air Blank | 0.000 | 08:56:09 |
| Subject Test | 0.183* | 08:59:17 |
| Air Blank | 0.000 | 08:59:49 |
| Air Blank | 0.000 | 09:00:16 |
| Subject Test | XXX** | 09:00:57 |
| Air Blank | 0.000 | 09:01:32 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006053
Core Version: 8105.48
01/30/2014 09:04:43

5 of 7 sc
TSH

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:05:11 |
| 0.100 Cal Check | 0.099 | 09:05:32 |
| Air Blank | 0.000 | 09:06:02 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
184
2.12.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/30/14 TIME 837
Intoxilyzer Serial # 80-006053 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006053
Core Version: 8105.48
01/30/2014 08:37:30

7 of 7 sc
TBH

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:37:59 |
| 0.100 Cal Check | 0.099 | 08:38:20 |
| Air Blank | 0.000 | 08:38:49 |

Preventative Maintenance
Performed -- Timer Reset