



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TJG  
6.3.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006053 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/22/14 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 4/24/14 and 5/22/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

5/29/14


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29th day of May, 2014.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

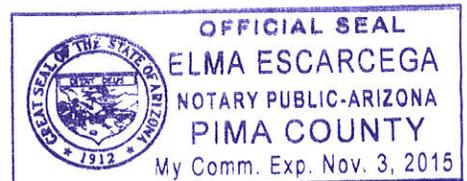

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS4  
6.3.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/22/14 TIME 1042

INTOXILYZER SERIAL # 80-006053 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
05/22/2014

10:46:29

2 of 7  
sc  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:47:29
Diagnostic Test	Pass	10:48:04
Air Blank	0.000	10:48:32
0.100 Cal Check	0.098	10:48:52
Air Blank	0.000	10:49:22
Subject Test	0.000	10:49:56
Air Blank	0.000	10:50:24
Wait	ABT*	
Air Blank	0.000	10:51:25

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
05/22/2014

10:53:35

3 of 7  
sc  
TSh

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:54:29
Diagnostic Test	Pass	10:55:04
Air Blank	0.000	10:55:32
0.100 Cal Check	0.099	10:55:52
Air Blank	0.000	10:56:22
Subject Test	RFI*	10:56:33
Air Blank	IPS**	10:56:45

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
05/22/2014

11:00:27

4 of 7  
SC  
Tsh

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:01:25
Diagnostic Test	Pass	11:02:00
Air Blank	0.000	11:02:27
0.100 Cal Check	0.098	11:02:48
Air Blank	0.000	11:03:18
Subject Test	0.256*	11:06:25
Air Blank	0.000	11:06:59
Air Blank	0.000	11:07:26
Subject Test	XXX**	11:07:55
Air Blank	0.000	11:08:30

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
05/22/2014 11:08:50

5 of 7<sub>sc</sub>  
Bh

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:09:19
0.100 Cal Check	0.098	11:09:40
Air Blank	0.000	11:10:09

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
sc  
TS4  
6.3.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/22/14 TIME 1042  
Intoxilyzer Serial # 80-006053 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
05/22/2014 10:42:53

*80 SC*  
*7 of 7 SC*  
*BSH*

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:43:23
0.100 Cal Check	0.099	10:43:44
Air Blank	0.000	10:44:13

Preventative Maintenance  
Performed -- Timer Reset